# S·E·Y·S·S·O

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This manual is intended to instruct you in the safe and correct operation of the device to avoid potential injury and damage.

**ALWAYS OBSERVE BASIC SAFETY RULES WHEN USING ELECTRICAL APPLIANCES, ESPECIALLY IN THE PRESENCE OF CHILDREN.**

**Dangers:**

1. Only use the USB cable and charging station included in the kit to charge your toothbrush.
2. Before connecting the charging station to the power supply, make sure that the voltage at the electrical socket is compatible with the specification of the device.
3. Do not use the toothbrush or charging station if it has been damaged in any way. If the charging station's power cord is damaged, stop using it immediately.
4. Do not attempt to repair, modify, or disassemble product components yourself. In case of damage, contact the customer service centre in your country or region (see the "Warranty and Service" section).
5. If the toothbrush will not be used for a long time, you want to clean it, or it is damaged, disconnect it from the charging station. If the toothbrush is not charging, unplug the charging station from the power supply.
6. If the toothbrush is fully charged, disconnect the charging station from the power supply.
7. Do not place metal objects on the charging station when charging.
8. Do not touch the USB cable or charging station with wet or damp hands.
9. Keep all parts of the product away from open flames and hot surfaces.
10. Product for indoor use only. Do not expose the product to external weather conditions.
11. When disconnecting the charging station from the electrical socket, pull the plug, not the cable.
12. Do not bend the USB cable and do not place it on sharp edges.

## Warnings:

1. Keep the charging station and USB cable away from water.
2. Do not immerse the charging station in water or any other liquid.
3. Make sure the charging station is completely dry before connecting it to the power supply.
4. To reduce the risk of electric shock:
  - Always disconnect the device from the power source after use.
  - Do not use while bathing.
  - For more detailed instructions on storing and using the device, as well as the included accessories, refer to the manufacturer (manufacturer's details are on the last page of the manual).
5. Do not immerse the toothbrush in water more than 1 metre deep or under running water for more than 30 minutes.
6. The toothbrush should be charged at a temperature between 0 and 40°C.
7. Use the product as intended. Follow the general safety precautions and battery recommendations in this manual. Any improper use may result in electric shock, burns, fire, and other hazards or injury.
8. Do not place any part of the product in a dishwasher, microwave oven, or on an induction cooker.
9. Do not expose any part of the product to direct sunlight or high temperatures (e.g. in hot cars or near hot cookers). Batteries may explode if overheated.
10. If the product shows signs of overheating, emits an unusual odour, or changes colour, or if charging takes significantly longer than usual, stop using the product, stop charging it, and contact customer service.
11. The toothbrush is equipped with a non-replaceable battery. Do not open the toothbrush to replace the battery.
12. If the battery is damaged, avoid contact with skin and eyes. In case of skin or eye contact, flush immediately with water and seek medical attention.

## Note:


1. The product may be used by children and persons with reduced physical, sensory, or mental capabilities only if they are properly supervised or have been instructed in the safe operation of the device and understand the risks associated with improper use. Cleaning and maintenance cannot be carried out by children without adult supervision.
2. Never use the product if you have concentration problems and feel drowsy.
3. Young children should not play with the product. They may swallow small parts, which can lead to choking.
4. The product is intended only for cleaning teeth and gums. It should not be used for any other purpose. If you feel pain or discomfort after using the toothbrush, or have any concerns about the effects of the device on your health, please consult your doctor.
5. If you have recently undergone surgery around the mouth

or gums, consult your doctor before using this product.

6. The toothbrush may cause your gums to bleed slightly in the first few days of use. If the bleeding does not stop within a week or is unusually heavy, or brushing is uncomfortable, contact your doctor.
7. Do not press the toothbrush with excessive force against the teeth.
8. Thoroughly clean the toothbrush after each use and make sure that there is no toothpaste remains on the handle or in the area of contact with the brush head.
9. For hygienic reasons, it is recommended that one exchangeable brush head should only be used by one person.
10. Replace the brush head once every 3 months, or more often if it wears out.
11. We recommend using only the replacement brush heads indicated by the toothbrush manufacturer.
12. The product complies with all standards for the operation of electromagnetic devices. If you have questions about the use of the toothbrush by people with artificial pacemakers or other implants, contact your doctor or the manufacturer of these devices.

## 2

### DESCRIPTION OF THE KIT COMPONENTS
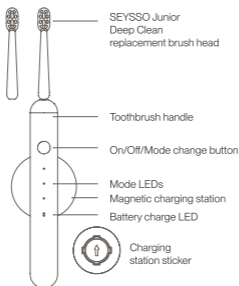
## 3

## STEP-BY-STEP INSTRUCTIONS FOR USE

**How to install and remove the brush head:**

**To install:** Follow the instructions in the picture. Make sure the head tip is aligned with the button and LEDs on the handle. Align them and place the head on the body. Make sure that the brush head is properly installed.

**To remove:** Follow the instructions in the picture. Grasp the brush handle with one hand and the head with the other. Pull vertically. We recommend replacing the brush heads every 3 months.

**Head compatibility**

SEYSSO Junior sonic toothbrush is compatible with SEYSSO Junior Deep Clean. We do not guarantee that the toothbrush will work properly with any other brush heads.

**How to install the charging station:**

Before installing the charging station, clean the surface intended for the sticker. It should be dry and degreased. Peel off the protective film from the sticker at the bottom of the station. Make sure to put the sticker properly on the surface. The arrow should be pointing upwards. Press gently on the sticking area. If the sticker holds tight, remove the stand to align the sticker and turn it clockwise to adjust the position.

Peel of the protective film from the sticker

Put the sticker on the surface - the arrow should be pointing upwards

Press gently on the sticking area

Turn the stand clockwise to adjust the position

**How to remove the charging station:**

Turn the charging station counterclockwise. When the charging station becomes loose, remove it by pulling it vertically.

**Note:**

To avoid damage, do not stick the charging station on whitewashed or wooden walls, wallpaper, or other uneven surfaces from which it can easily peel off. It should be stuck on flat and smooth surfaces, such as glass or ceramic tiles.

**Note:**


The charging base has a type-C charging port. The charging station should be installed in a dry environment. Ensure that the station port is dry when charging.

**Starting and stopping**

Press the button briefly to start the toothbrush. Press the button again after at least 5 seconds of operation in the current mode to switch the toothbrush off.

**Modes**


Press the button within 5 seconds of starting or changing the mode to change the mode.

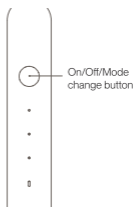**Memory**

The current mode is stored after 5 seconds of use. The next time you start the toothbrush, it will set the mode you selected previously.

**Modes:****How to set the mode**

Press the button briefly to change the mode. LED on the toothbrush handle will indicate the selected mode.

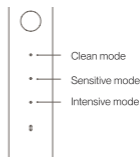**Clean**

Rapid removal of plaque and stains with 76,000 sonic movements per minute. Recommended for daily use.

**Sensitive**

Low-frequency vibrations – 56,000 sonic movements per minute. Suitable for people with sensitive gums or teeth. Also suitable for people who use a sonic brush for the first time.

**Intensive**

High-frequency vibrations – 96,000 sonic movements per minute. Allows you to thoroughly remove stains and helps polish teeth.

## 4

RECOMMENDATIONS  
ON BRUSHING TECHNIQUES

1. Put the toothpaste on the brush head of the toothbrush.
2. Divide your mouth into four sections: upper right, upper left, lower right, and lower left. Brush each section for 30 seconds. The toothbrush will signal when to move to another section.
3. Place the brush head in your mouth and press the button to switch the brush on.
4. Keep pressing the button to switch between different modes: Clean (for cleaning), Sensitive (for sensitive teeth) and Intensive (for removing stains).
5. Clean all tooth surfaces one at a time: cheek-sides, chewing surfaces, and lingual surfaces. Use slow movements, checking how firmly you press the brush against your teeth.
6. The chewing surfaces should be cleaned at a 90 degree angle. The facial and lingual surfaces should be cleaned at a 45 degree angle. If you press the toothbrush against the teeth too hard, the head movement may become weaker during cleaning.
7. Rinse the brush head and the handle under running water. Dry them and place them in a dry place.

**Note:**

To avoid splashing liquids, do not take the toothbrush out of your mouth while it is running.

**Replacement brush heads**

After each use of the toothbrush, the brush head should be removed from the handle, rinsed with running water, and allowed to dry.

**Handle**

Rinse the toothbrush handle under running water to remove any remaining toothpaste and wipe dry with a soft towel.

**Charging station**

Unplug the charging station and clean with a damp towel. Allow to dry.

**Note:**

Do not use water above 40°C and acidic or strongly alkaline cleaning fluids to clean the kit components.

**Inductive charging**

Your new toothbrush is partially charged. The toothbrush needs to be recharged when the battery is empty. If the charging indicator flashes red, this means that the toothbrush needs to be recharged. While charging, the battery charge LED will flash red. When fully charged, the LED will turn white. For charging, use a USB charger (adapter/plug) with a standard voltage of 5V. A charger with a standard charging current of 1A is recommended.

**Note:**

To prolong the life of the battery, charge your toothbrush only when prompted to do so (when the charging indicator LED is flashing red). Keep charging until the battery charge LED turns white.

**Disposal**

The device can only be disassembled when it is intended for disposal. Otherwise, you may destroy the product. The product has a built-in rechargeable battery. Do not dispose of the battery with municipal waste. Please dispose of the product at an official electrical waste collection centre.

## 6

## TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The vibration of the brush seems weaker.	The replacement head may have been incorrectly applied to the handle.	Make sure the head is applied correctly.
		If the problem persists, replace the head with a new one.
The battery life is very short when charged.	The charging process may be ineffective due to incorrect positioning of the toothbrush on the charger.	Make sure that the toothbrush is placed upright while charging.
		If the problem persists, the likely cause is a spent battery. Contact Customer Service.
The charging indicator is not lit when the toothbrush is connected to the charging station.	If the toothbrush has not been used for more than six months, it enters sleep mode.	After 6 to 20 minutes, the toothbrush should awaken from its sleep state and start charging the battery. This is normal.
		If the toothbrush still does not charge, contact Customer Service.
The toothbrush does not start.	Check if the power supply is sufficient.	Ensure that the toothbrush has been properly charged.
		If the toothbrush still does not start after charging, contact Customer Service.

## 7

SEYSSO JUNIOR SONIC  
TOOTHBRUSH SPECIFICATION

PRODUCT NAME	SEYSSO JUNIOR SONIC TOOTHBRUSH
Voltage	5V/1A
Power	3W
On/Off Switch	Button
Display	LED
Displayed Information	Mode + Low battery indication
Usage Time	120 - 130 minutes
Charging	Inductive charging
Charging Time	about 10 hours
Charging current	80±15mA
Battery	800mAh
Discharged Battery Indicator	Flashing diode
Waterproof Rating	IPX7

**Important: Use of the appropriate power supply**

Please ensure that the power supply used has output parameters compatible with those specified on the label of your toothbrush. Using a power supply with incompatible parameters may lead to improper functioning of the device and is absolutely unacceptable. Failure to adhere to this rule may result in damage to the device, as well as loss of warranty. Take care of your device by using only recommended accessories and power supplies.

The warranty is valid in Poland and in the territories of official distributors. The warranty covers damage related to manufacturing defects or improper workmanship of the product. Keep the proof of purchase for the duration of the warranty. During the warranty period, you can file a claim if the product was used correctly and for its intended purpose.

#### THE WARRANTY DOES NOT COVER:

1. Surface damages resulting from normal use.
2. Damages resulting from:
  - a) accidental dropping or impact, improper use, abuse or commercial use, neglect, device modification,
  - b) improper maintenance, neglect of the product,
  - c) connection to an incorrect power source,
  - d) repairs or attempted repairs by an unauthorized party by the manufacturer.

Any attempts to open or dismantle the device, its parts, or accessories will void the warranty. Warranty claims must be supported by proof. The claim must be made within the warranty period. Scratches or abrasions on the surface of the handle will be treated as superficial damage resulting from normal use. As such, they will not be covered by the warranty. The warranty applies to components that affect the functioning of the toothbrush. Proper care of your sonic toothbrush significantly increases its lifespan. Failure to follow recommendations may affect the loss of warranty.

INSTRUCTION FOR MODELS: SK.ST.010.GRN.1, SK.ST.010.PNK.2.

**PRODUCER:**

WM BRANDS SP. Z O.O.

UL. DZIUPIŁI 12

02-454 WARSZAWA, POLAND

MADE IN P.R.C.

**CONTACT:**

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WEBSITE: WWW.SEYSSO.PL

Water Resistance (Toothbrush Handle): IPX7. Protection against the effects of immersion in water up to a depth of 1 meter for 30 minutes. The toothbrushes are equipped with a built-in lithium-ion battery.

This product meets the requirements of the so-called "New Approach" Directives of the European Union (EU).

The packaging is made of cardboard. It is recyclable.

The instruction manual is made of paper. It is recyclable.

This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU with the symbol of a crossed-out waste bin. This marking indicates that the equipment, after its period of use, should not be disposed of with other household waste. The user is obliged to hand it over to entities collecting waste electrical and electronic equipment, including local collection points, stores, and municipal units, which create an appropriate system allowing for the disposal of this equipment. Proper handling of waste electrical and electronic equipment helps avoid harmful consequences for human health and the environment due to the presence of hazardous components and improper storage and processing of such equipment.

The device is marked with the symbol for Class III protection for electrical devices.

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