S·E·Y·S·S·O

OXYGEN O-SONIC
TOOTHBRUSH

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EN / PL

The instruction applies to model with
the following number PSTO01

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Before using this product, please
read carefully this user's manual and
keep it properly for later reference.

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
1 Attention:

Please read the manual carefully before using.

Warning:

1. Before every use check if device is damaged or worn. Do not use if damaged or worn.
2. Do not charge it if the charger has any damage.
3. Stop using the device immediately if it has any signs of damage.
4. The charger cannot be used outside or nearby the extremely hot surfaces.
5. Do not use while you are taking a bath or shower.
6. Please consult your dentist before you use if you have had undergone an operation of your mouth, gums or teeth over the past 2 months.
7. Please consult your dentist if (oral) bleeding occurs especially when it lasts over 1 week after use.
8. Use only original brush heads.
9. Please consult your doctor or the manufacturer regarding use while having implanted a pacemaker or other similar devices.
10. The product is only suitable for cleaning teeth, gums and tongue. Do not use it for other purposes. Please stop using the device and contact the doctor if you feel unwell or painful after use. Please consult the doctor before (first) use if you have any concerns about your health.
11. The product is improper for children, babies, and operators with reduced physical, sensory or mental capabilities unless they are supervised or guided by capable person.
12. Store the product beyond children's and babies' reach. Do not use as toy.
13. The power cord should not be replaced. If the cord is damaged the appliance should not be used.

2 Specifications:

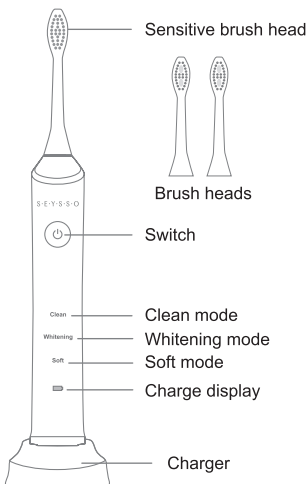Product Name:	Sonic Toothbrush
Working Voltage:	100-240V/AC
Rated Power:	3W
Switching Pattern:	Touch Switch
Display:	LED
Display Content:	Clean/Whitening/Soft + Battery Indicator
Vibration Frequency:	62000-96000 times per minute
Working Time:	About 250 minutes
Charging Mode:	Inductive Charging
Charging Time:	About 15 hours
Charging:	Current 55mA
Battery Specification:	Lithium Battery/14500/600mAh
Power Warning:	When the power is low, the power indicator light flashes
Waterproof Grade:	IPX 7

3 Troubleshooting:

Issue	Possible Reasons	Solution
The brush head movement is weakening.	(1) The toothbrush is pressed too hard against the teeth. (2) The brush head is attached too close to handle.	(1) Decrease the pressure. (2) Take the brush head off. Place it again on the pin. Be aware the gap between handle and brush head is needed.
The toothbrush works for a very short period of time even if it has been charged.	(1) The handle was not put on the charger properly. (2) The battery is no longer charged.	(1) Place the handle properly, continue charging. (2) Charge the toothbrush before using if you have not used it for over 6 months. (3) Contact the After-Sale Service.

4 Charging:

1. The battery comes partially charged. Keep using until complete discharge. Then use charger.
2. If you have not used the toothbrush for over 6 months charge it for over 15 hours. Fully charged it would be able to run for up to 60 days without recharging.
3. To prolong the battery life charge only if the battery has no charge left. Lithium-Ion battery life can be shortened if done otherwise.
4. If the charge indicating light blinks for several seconds continuously and the device stops working the battery has low power.
5. Put the handle on the charger, confirm that it is connected to the power supply. The indicating light should start to blink green showing that the charging is in progress.
6. The indicating light keeps blinking as the battery charges. When battery is fully charged, the indicating light should lit constantly.
7. When left on the charger after 20 minutes of being fully charged the toothbrush will enter recharging mode automatically.
8. Fully discharge the toothbrush at least once a year. Unplug the charging stand and exhaust the battery by using the toothbrush until it has no power left.

5 Maintenance:

Clean the brush head:

1. After each use clean the bristles of the brush head and let them dry.
2. Do not pull (out) the bristles.

Clean the handle:

1. Clean it with water let it dry.
2. Do not immerse it in water.
3. Do not use hot water (over 40 degrees Celcius) or acidic cleansing solution.

Clean the charger:

1. Unplug the charger before cleaning.
Cleaning plugged device may cause electric shock.
2. Clean charger with soft cloth.
3. Clean the plug with dry cloth in order to avoid the short circuit.

Change the brush head:

Dentists advise brush head change 3 months.


6 Proper brushing time:

1. Dentists advise that each time you should brush your teeth for 2 minutes.
2. When using, you will feel pause every 30 seconds, which reminds you to clean the four areas in your mouth evenly in order to achieve the best cleaning effects.


7 Using, step by step:

Step 1: For your first use please charge the toothbrush for 15 hours.


Charge for 15 hours

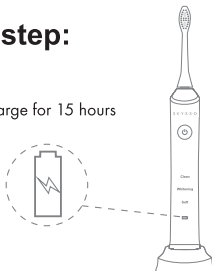Step 2: Attach the brush head onto the pin carefully. Pay attention that the gap must be left between the brush head and the handle.

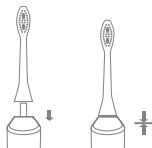Step 3: Squeeze moderate amount of toothpaste on the bristles, soak it, put it into your mouth, and start the toothbrush.

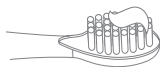Step 4: Align the brush head along the teeth and gum line at an angle of 45.

Step 5: Hold the toothbrush gently and move it slowly against the gum line. Do not press.

Step 6: Rinse your mouth with clean water and wash the brush head afterwards.


8 3 professional brushing modes:

The product has 3 brushing modes: clean, whitening and soft to customize the brushing experience especially for you.

You can switch the modes in order to care and beautify your teeth in an all-round manner.

- Clean
- Whitening
- Soft

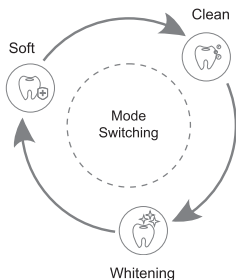Brushing Mode Memory:

1. Press the on/off button and the toothbrush begins to work. If any brushing mode works for over "5" seconds it will be memorized when on/off switch is pressed.
2. Press the on/off button and the toothbrush begins to work. If you shut it down by switching the modes, when restarted - it will return to the default "Clean" mode.

9 Notice:

1. Do not share the brush head with your other users.
2. Do not intertwine the power cord line with the adapter tightly when you store the charger.
3. Do not put metal objects (such as coins, clips and so on) on the charger.
4. To avoid fire unplug the charger when it won't be used in foreseeable future.

10 Notes:

It is possible to have irritated gums or slight bleeding when using the toothbrush for first few times. This is because the gingiva is not used to such stimulation. These symptoms will usually cease within 1 or 2 weeks. If they last over 2 weeks that might indicate some health issues. Under this circumstances, you should stop using the toothbrush and consult the dentist.

11 Warning:

1. Do not immerse the charger in water or any other liquids.
2. After cleaning ensure that the charger plug is completely dry before plugging it back to the power supply.
3. This product contains built-in rechargeable battery. Do not put it into fire or heat it. Do not charge, use or leave it exposed on high-temperature.
4. Do not strike or flop this product. Do not use if damaged by such activity.
5. Do not remould, disassemble, or fix this product by yourself.
6. When charging keep the product away from the strong electromagnetic field (like the induction cooker).

12 Replace brush head:

SEYSSO warranty will be invalid if it is found that the handle is damaged because non-SEYSSO replacement brush heads are used.

It is not recommended to use non-SEYSSO replacement brush heads, because SEYSSO cannot control the quality of other replacement brush heads.

SEYSSO cannot guarantee that non-SEYSSO replacement brush heads fit the handle.

All SEYSSO replacement brush heads are printed with the logo of SEYSSO and meet the high-quality standard of SEYSSO.

13 Warranty and services:

You will enjoy two-year quality warranty for this product since the day of purchasing. In the warranty period, we will be responsible for repairing or replacing the complete device if malfunction is caused by any design defect, but the warranty does not cover the damages caused by inappropriate use of the product.

Quality warranty – restrictive clause

Quality warranty is only applicable to the products manufactured by SEYSSO.

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Explanation:

Please carefully fill in the warranty card and keep it appropriately for later service. You can enjoy two-year free warranty services for the product since the day of purchasing. You cannot enjoy free services if any of the following circumstances happens:

1. The product is damaged because of inappropriate use or improper storage;
2. The product has been disassembled or repaired by others rather than SEYSSO special repair centers;
3. You do not have effective proof of purchase;
4. Damages are caused by instable voltage or voltage beyond the specialized range;
5. The product is not in the warranty period or beyond the coverage of warranty.

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Kontakt: -----

Fault description:

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The instructions apply to models with the numbers: PSTO01

MANUFACTURER:
WM BRANDS SP. Z O.O.
UL. DZIUPI 12 02-454
WARSAW, POLAND
Made in P.R.C.

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Power: 3W

Water resistance (handle): IPX7

Protection against the effects of immersion in water at a depth of up to 1 meter for 30 minutes
Toothbrushes have a built-in lithium-ion battery.

The product meets the requirements of the European Union (EU) "New Approach" directives.

The packaging and instruction manual are recyclable.

This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU with the symbol of a crossed-out waste bin. Such marking informs that this equipment, after its use, cannot be disposed of together with other household waste. The user is obliged to hand it over to entities collecting waste electrical and electronic equipment. Collectors, including local collection points, stores, and municipal units, create an appropriate system allowing for the disposal of this equipment. Proper handling of waste electrical and electronic equipment helps to avoid harmful consequences for human health and the natural environment, resulting from the presence of hazardous components and improper storage and processing of such equipment.

Symbol of Class II protection designation for electrical equipment.

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