S·E·Y·S·S·O

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
The instruction applies to models with the following numbers /

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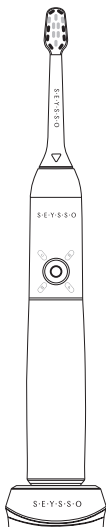EN SONIC TOOTHBRUSH

DETAILED INSTRUCTION MANUAL
BEFORE USE, PLEASE READ THE INSTRUCTION MANUAL.
KEEP THE MANUAL FOR FUTURE REFERENCE

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READ BEFORE USE!

KEEP THIS MANUAL


- For domestic use only. 
- Place the toothbrush on the charger and connect the charger plug to an outlet (socket specification: 100-240V~50-60Hz).
- The toothbrush is equipped with a waterproof handle. It is safe and designed for use in the bathroom. The toothbrush has an IPX7 protection rating, meaning it is resistant to accidental immersion in water up to a depth of 1 meter, for up to 30 minutes.
- Remember: the toothbrush is not designed to operate underwater!
- The product is not intended for self-repair.
- Store the charger in a dry and safe place.
- Straighten the cable before use.
- In case of failure, do not disassemble the device or replace any components as there is a risk of electric shock. Contact SEYSSO.
- A damaged power cord must be replaced with a new one. Contact SEYSSO.
- If dropped, replace the head before next use. Even if the damage is not visible.

It is not recommended to use the device by children under the age of 8. Children and persons with limited sensory, motor, or mental capacities should use the device under the supervision of a caregiver or should be properly trained on the correct and safe use of the device. Children should not play with the device. Unsupervised children should not perform cleaning and maintenance of the device.

WARNING ABOUT DANGER

- To reduce the risk of electric shock:
- always disconnect the device from the power source after use.
- do not use while bathing.
- for more detailed instructions regarding storage and use of the device, as well as included accessories, refer to the manufacturer.
- To reduce the risk of ignition, short circuit, or injury:
- disconnect the charger from power after charging is complete.
- the product is not a toy. Exercise special caution if the toothbrush is used by or near children.
- use according to its purpose. Only use accessories recommended by the manufacturer.
- do not use the charger if the power cord or plug is damaged. If the toothbrush does not work properly after being dropped, left outdoors, or immersed in water, return it to the place of purchase.
- do not stretch the charger's power cord, trap it in doors, wrap it around sharp edges. Do not turn on the toothbrush while charging. Keep device components away from heated surfaces.
- do not disconnect the charger from the power by pulling on the cable. Use the plug.
- do not handle the plug with wet hands.
- never use the product if you have problems concentrating or feel drowsy.
- do not use the device outdoors or in places where it can be exposed to aerosols or pure oxygen.
- connect the device to an outlet. Familiarize yourself with its specification.

2 Construction of the toothbrush

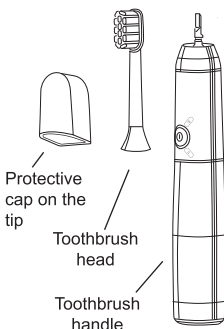3 Charging the toothbrush before first use


1. Before the first charge, the toothbrush should be completely discharged. Then charge the toothbrush to full. The first charge should last 32 hours. During this time, the toothbrush can be removed from the charger for the duration of brushing teeth.
2. The LED ring on the toothbrush flashes in four different colors during charging. The LED ring turns off 2 minutes after the toothbrush is fully charged.
3. Each subsequent charge should last 32 hours or until the LED ring goes off.

NOTE: Before charging, the toothbrush should be turned off. The charger base should be placed on a flat and stable surface, and the hole at the bottom part of the toothbrush should be aligned with the locating pin of the charger base. The toothbrush should be placed on the charger, and the plug connected to a power outlet.

NOTE: Before starting charging, check if the output current specification is compatible with 100-240 V AC 50 Hz/60 Hz 1.4 W. Otherwise, do not connect the charger to the power supply.


4 Instructions for using the toothbrush

Step 1: Attach the replaceable head to the toothbrush handle

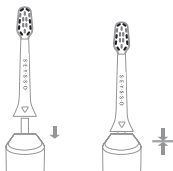• The diagram shows the correct way to attach the head to the toothbrush.

Step 2: Then apply the appropriate amount of toothpaste to the bristles. We recommend using toothpastes dedicated for sonic toothbrushes (Sonic Friendly). We recommend SEYSSO toothpastes, which thanks to their ideal composition and consistency fully exploit the potential of sonic toothbrushes.

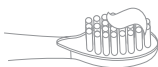Step 3: Before turning on the toothbrush, place the head in the mouth. Position the toothbrush at the teeth, so that its bristles are set at a 45-degree angle to the gum line and turn on the toothbrush by pressing the function button for 3 seconds.

Step 4: Briefly press the function button during the toothbrush's operation to select the appropriate tooth cleaning program. Each press of the button switches the device alternatively between four programs and changes the color of the LED ring.

Mode	LED ring color	Symbol on the toothbrush	Function
Clean	Yellow	C	Tooth cleaning
White	Blue	W	Whitening
Polish	Red	P	Additional polishing of front teeth
Massage	Green	M	Gum massage and tooth cleaning

Step 5: Slowly move the toothbrush along the gum line, without pressing too hard on the teeth. The toothbrush has a built-in automatic, 2-minute timer, which helps adhere to the dentist-recommended tooth brushing time. Thanks to 30-second intervals, indicated by a brief pause in operation, you can evenly divide the tooth brushing time across all quadrants of your teeth.

Step 6: The device will turn off automatically after 2 minutes. To turn off the toothbrush during operation, hold the button for 3 seconds.

Step 7: After brushing your teeth, remove the head from the toothbrush handle, rinse with water and dry before storing. Remember to remove the head after each use of the toothbrush. When the toothbrush battery is depleted, place the toothbrush on the charger base and charge according to the instructions described above.

Step 8: When the battery charge level falls below 3.3V, the LED ring will start flashing slowly. When it falls below 3V, the LED ring will start flashing quickly and the toothbrush will automatically turn off.

5 Replaceable heads:

The set includes a SEYSSO Planet XCleaning head. Use only original SEYSSO heads dedicated for this toothbrush. Any damage caused by using non-original heads will not be covered by the manufacturer's warranty.

6 Cleaning and maintenance

- After brushing your teeth, rinse the toothbrush under running water.
- Remove the head from the handle and rinse it under running water.
- Wipe the handle before placing it on the charger.
- Before cleaning the charger, make sure it is unplugged. Use only a dry cloth to wipe the charger surface!

7 Two-year warranty

The manufacturer provides a warranty for two years from the date of purchase. The warranty is valid in Poland and in the territories of official distributors. The warranty covers damages related to manufacturing defects or improper manufacture of the product. Keep the proof of purchase for the duration of the warranty. During the warranty period, you can submit a claim if the product was used properly and according to its intended use.

The warranty does not cover:

1 Surface damage resulting from normal use.

2 Damages resulting from:

- a) accidental dropping or impact, improper use, misuse or commercial use, neglect, modification of the device,
- b) improper maintenance, carelessness about the product,
- c) connection to an improper power source,
- d) repairs or attempted repairs made by an unauthorized party by the manufacturer.

Any attempts to open or dismantle the device, its parts, or

accessories will void the warranty. Warranty claims must be supported by evidence. The claim must occur within the warranty period. The material covering the handle is delicate, so care must be taken during use. Scratches or wear on the handle surface will be treated as surface damage resulting from normal use. Thus, they will not be covered by the warranty. The warranty applies to components, which affect the functioning of the toothbrush. Proper care of your sonic toothbrush significantly extends its lifespan. Non-compliance with recommendations can affect the loss of warranty.

This manual applies to models numbered: SPST10, SPST20

MANUFACTURER: WM Brands sp. z o. o. ul. Dziupli 12 02-454
Warsaw, Poland
Made in P.R.C.

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Website: www.seysso.pl

Charger input voltage: 100-240 V~, 50/60 Hz

Input current: 35mA Power: 1.4W

Water resistance (handle and charger of the toothbrush): IPX7.

Protection against effects of immersion in water up to 1 meter deep for 30 minutes. The toothbrushes have a built-in lithium-ion battery.

 The product meets the requirements of the so-called "New Approach" Directives of the European Union (EU).

 The packaging is made of cardboard. It is recyclable.

 The manual is made of paper. It is recyclable.

 This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and Directive WEEE 2012/19/EU with the symbol of a crossed-out waste bin. Such marking informs that this equipment, after its period of use, cannot be disposed of together with other household waste. The user is obligated to hand it over to entities managing the collection of waste electrical and electronic equipment. Those managing collection, including local collection points, stores, and municipal units, create an appropriate system enabling the disposal of this equipment. Proper handling of waste electrical and electronic equipment helps avoid harmful consequences for human health and the natural environment resulting from the presence of hazardous components and improper storage and processing of such equipment.

 Symbol indicating Class II protection for electrical devices.

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