
Feelo.

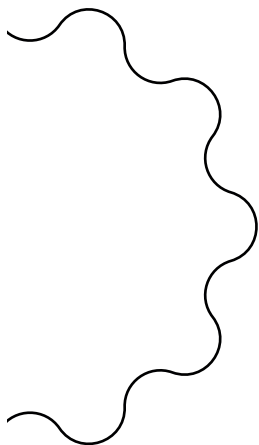PRO

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SET ELEMENTS

- 1 brush handle
- 1 replacement head
- 1 travel case
- 1 charging cradle with integrated USB cable (without adapter)

TOOTHBRUSH FEATURES AND FUNCTIONS

- Wide speed range: 48,000-66,000 sonic movements per minute ensure effective cleaning of teeth adjusted to the user's needs.
- 4 work programmes: Clean (daily cleaning), White (white), Polish (polish), Gum Care (gum care). 3 brushing intensities: high, medium and low.
- Timer – 2-minute automatic timer, reminding every 30 seconds to move on to the next quadrant of the mouth.
- Powerful lithium-ion battery – 1 charging is enough for about 60-70 days of use at low brushing intensity.
- Headbrush replacement indicator – reminds you when to replace the head with a new one.
- IPX8 high water-resistance rating. This means that the toothbrush is protected against the effects of longer submersion in water.
- The replaceable brush head has thick bristles, which ensures thorough cleaning of teeth and hard-to-reach places in the mouth. The rear rubber part of the head is used to clean the tongue.

SELECTION OF THE TOOTHBRUSH FUNCTIONS

Switching on the toothbrush	Briefly press the 'Smart button' to switch on the toothbrush.
Switching off the toothbrush	Briefly press the 'Smart button' after 3 seconds from switching on the toothbrush to switch it off.
Changing the work programme	To change the work programme, briefly press the 'Smart button' within 3 seconds after switching on the toothbrush.

SELECTING BRUSHING INTENSITY (STRENGTH)

To select your preferred brushing intensity, press and hold longer the 'Smart button' after 3 seconds from switching on the toothbrush. The 'Strength' indicator will glow with different color depending on the brushing intensity selected.

- **Low cleaning intensity:** 'Strength' indicator – glows green. Strength of the cleaning appropriate for sensitive teeth and gums.
- **Medium cleaning intensity:** 'Strength' indicator – glows blue. Strength of the cleaning appropriate for daily tooth brushing.
- **High cleaning intensity:** 'Strength' indicator – glows pink. Strength of the cleaning appropriate for removing persistent discolorations.

RESETTING THE HEAD REPLACEMENT REMINDER

Press and hold the 'Smart button' for 5 seconds when the toothbrush is switched off.

(When the toothbrush is switched on and off 180 times, the head replacement indicator will start flashing. Replace the brush head with a new one and reset the head replacement reminder.)

RESTORING THE FACTORY SETTINGS

Pressing the 'Smart button' for 5 seconds while charging the appliance restores the factory settings. After restoring the factory settings, the head replacement indicator will flash 3 times. (The toothbrush memory function remembers the last selected work programme and cleaning intensity, therefore it is not necessary to set them again each time after switching on the toothbrush. The settings will be retained until you change them.)

DESCRIPTION OF THE TOOTHBRUSH WORK PROGRAMMES

PROGRAMMES	
CLEAN	Programme for daily tooth brushing. High brushing strokes frequency ensures effective removal of dental plaque.
WHITE	The variable brushing strokes frequency helps remove discolorations, restoring natural white teeth.
POLISH	The high brushing strokes frequency with a small deflection polishes and brightens the teeth, making them shiny and smooth.
GUM CARE	Softer strokes of the toothbrush ensure gentle cleaning of the teeth and strengthening massage of the gums, improving their condition.

DESCRIPTION OF THE TOOTHBRUSH FUNCTIONS

FUNKCJE	
STRENGTH	By pressing the button you can switch between 3 cleaning intensities.

	When the power is low, the charging indicator will start flashing and the toothbrush handle will vibrate 3 times to remind you to charge the toothbrush.
	The head replacement indicator flashes, reminding you to replace the cleaning head with a new one. It is recommended to replace the head regularly every 3 months.
3 sec	For the first 3 seconds after switching on, the toothbrush vibrates more gently to prevent toothpaste splattering.

POSSIBLE COMBINATIONS OF WORK PROGRAMMES

The Feelo Pro toothbrush has four work programmes and three brushing intensity levels which can be freely combined to create your own personalised toothbrushing mode, adjusted to individual needs of the mouth. The toothbrush offers 12 different cleaning combinations.




The picture below shows illustratively possible combinations of the brushing programme and cleaning intensity, where:

- **RED** means high brushing intensity,
- **BLUE** means medium brushing intensity,
- **GREEN** means low brushing intensity.

Each cleaning programme lasts 2 minutes.

USING THE TOOTHBRUSH

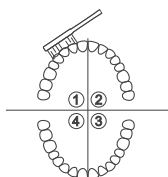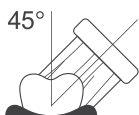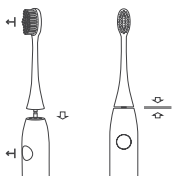- Fitting the brush head: place the front of the brush head in the same direction as the 'Smart button' and press the brush head vertically against the handle. There should be a small gap between the brush head and the handle.
- Wet the head tip with water and apply an appropriate amount of toothpaste.
- Place the brush head so that its bristles are at a 45° angle to the gum line. Make sure that the bristles also reach the gums. After switching on the toothbrush, select the work programme and its intensity. Remember to thoroughly clean each surface of your teeth (the inside, outside and chewing surfaces).
- Brush each quadrant of your teeth for 30 seconds. The timer divides the tooth brushing time into four 30-second intervals. A slight change in the vibration signals when it is time to move on to cleaning the next area of your mouth.
- The toothbrush switches off after 2 minutes of work, i.e. after the time recommended by dentists.

- Remember to brush your teeth gently. Brushing too hard can damage your gums and cause the abrasion of the tooth enamel.
- To avoid toothpaste splattering, always switch on the toothbrush only after you have placed it against your teeth, and only remove it from your mouth after brushing your teeth.
- With each subsequent use, the toothbrush will switch on on the previously set programme and intensity.

CHARGING THE TOOTHBRUSH

- The first charging of the toothbrush and each subsequent charging lasts 14 hours.
- To charge the toothbrush, place it on the charging base. Connect the charging station to a power source using a USB charger. The charging indicator will start flashing.
- The operating time of the toothbrush on a single charge is about 60-70 days at low brushing intensity.
- The toothbrush should be charged if it has not been used for more than 3 months.
- To prolong the service life of the battery, charge your toothbrush only when it indicates such a need. Continuous charging or storing a fully discharged toothbrush may reduce the battery performance.
- When you want to charge the toothbrush, remember not to place a wet handle on the charger.

LOW BATTERY LEVEL 20%	The charging indicator blinks and the toothbrush handle vibrates 3 times
CHARGING	The charging indicator flashes
FULLY CHARGED BATTERY 100%	The charging indicator glows continuously

SOLVING PROBLEMS

PROBLEM	SOLUTION
The toothbrush won't turn on.	Charge your toothbrush.
	If the toothbrush still does not start after the charging, contact the customer service department.
The vibration of the brush decreases.	Charge your toothbrush or replace the brush head. The brush vibrates differently depending on the selected work program.
The charging light does not light up when the toothbrush is connected to charge.	If the toothbrush has not been used for a long time, it enters a sleep mode. After a short time the toothbrush should awaken from the sleep mode and begin charging the battery. This is a standard behaviour.
	If the toothbrush still does not charge, contact the customer service department.
Tickling when brushing teeth.	The movements of the sonic toothbrush can cause a tickling sensation, especially the first time you use it. This impression will disappear with time.

SPECIFICATION

PRODUCT NAME	Feelo Pro Sonic Toothbrush
Rated voltage	5 V
Power	2.6 W
Switch type	Button
Display	None
Information displayed	Work programmes, brushing intensity, battery charging, head replacement indicator.
Operation time	Up to 60-70 days at low brushing intensity.
Way of charging	Wireless inductive charging
Charging time	About 14 h
Battery	Li-Ion 2000 mAh
Battery discharge signalling	Flashing of the charging indicator and vibrations
Water resistance class	IPX8

CLEANING AND MAINTENANCE

- After each brushing, remove the head from the toothbrush. Then rinse the handle and the head thoroughly under running water.
- Also clean the base of the handle to remove the toothpaste sediment and residue.
- Make sure that the handle put away is not wet. Wipe the handle dry. Store the toothbrush in a dry place at room temperature.

WARNINGS

1. Do not disassemble, modify or attempt to repair the product yourself.
2. If you have undergone oral surgery or have been diagnosed with severe periodontal disease, consult your physician before use.
3. Children and people unable to operate the product on their own should only use it under appropriate adult supervision.
4. Do not immerse the product in water or other liquids for a long time.
5. Do not expose the product to high temperatures or sunlight for a long time.
6. If the power cable of the charger is damaged, replace the charger.
7. Do not press the toothbrush hard against the gums and teeth.
8. The toothbrush is intended only for cleaning the mouth, do not use it for other purposes.

TWO-YEAR GUARANTEE

The Feelo toothbrush is covered by a 2-year guarantee from the date of purchase. The product guarantee is provided by Shop Dent Wojciech Madej with its registered office at Dziupli Street 12 in Warsaw. The guarantee is valid in countries in which the product is sold by Shop Dent Wojciech Madej or an authorised distributor.

The guarantee covers damage related to manufacturing defects or incorrect workmanship of the product. Under the guarantee, the guarantor undertakes to repair or replace a damaged element. Keep the proof of purchase for the duration of the guarantee period. Within the guarantee period, the buyer can make a claim to the place of purchase only if the product has been used properly for its intended purpose.

The guarantee does not cover:

- 1. Surface damage resulting from normal use:**
Scratches or abrasions on the surface of the of the handle will be treated as surface damage resulting from normal use and will therefore not be covered by the warranty. The guarantee applies to the components that affect the functioning of the toothbrush.
- 2. Damage resulting from:**
 - Accidental dropping or impact, misuse, abuse or commercial use, negligence or modification of the appliance.
 - Improper maintenance, lack of care of the product.
 - Connecting to an incorrect power source.
 - Repairs or attempted repairs by a party not authorised by Feelo. Any attempts to open or dismantle the appliance, its parts or accessories will render the guarantee invalid. Guarantee claims must be supported by evidence. A claim must be made within the guarantee period.

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