
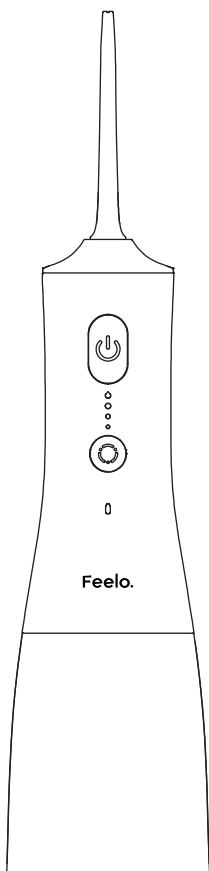
# Feelo.
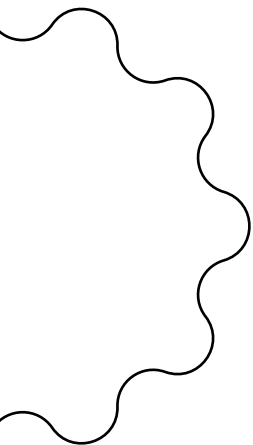
## PRO

IRYGATOR  
WATER FLOSSER  
ІРИГАТОР

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

## FEATURES AND FUNCTIONS OF THE IRRIGATOR

- Wide range of water pressure - from 2.07 to 7.58 bar.
- 4 operation programmes: Soft (gentle cleaning for those with sensitive gums), Standard (basic programme for daily irrigation), Strengthen (higher pressure programme for deep cleaning of hard-to-reach areas in the mouth), Pulse (for simultaneous cleaning of teeth and gum massage).
- Charging via the USB power cable USB-C ensures easy and convenient use of the irrigator at home or on the go.
- Powerful lithium-ion battery with a capacity of 2000 mAh - one full charge of the irrigator lasts for approx:
  - 90 days of use in mode Soft,
  - 80 days of use in mode Standard,
  - 60 days of use in mode Strengthen,
  - 70-80 days of use in mode Pulse.Assuming the irrigator is used by one person twice a day.

- Fast charging in 4 hours.
- The slim and ergonomic handle fits well in the hand for comfortable use.
- High degree of water resistance IPX7. This means that the irrigator is resistant to accidental immersion in up to 1 metre of water for up to 30 minutes.

## DESCRIPTION OF THE IRRIGATOR

1. Top cover
2. Tip lock release button
3. Flow opening
4. Anti-slip coating
5. Irrigator body/handle
6. USB-C type charging port and rubber cap protecting the port
7. Opening and closing the water container
8. Water container flap
9. Water container
10. On/Off Button
11. Change mode button
12. Charging indicator
13. Standard irrigator nozzle
14. Pulse mode
15. Strengthen mode
16. Standard mode
17. 1Soft mode
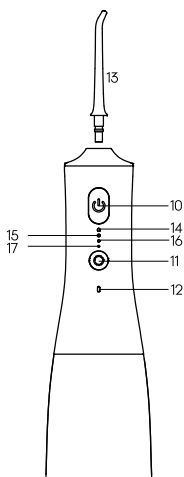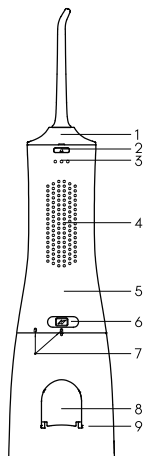
## NOTE

There may be residues of distilled water in the new unit (in the container, suction line, etc.), which was used for the performance test of the unit. It should not cause any concern when using the irrigator.


## CHARGING THE DEVICE (FIG. 1)

- Remove the cap covering the charging port and make sure there is no water residue inside the socket.
- Connect the irrigator to a power source using the cable USB-C to your computer or to the charger USB, to start charging.
- The charging indicator flashes when the device is charging. When the battery is fully charged, the charging indicator lights up white for 2 minutes and then switches off automatically.
- When charging is complete, put the cap back on the charging port.

## NOTE

If the charging indicator flashes during or after irrigation, this indicates that the battery needs charging. It takes approximately four hours to fully charge.
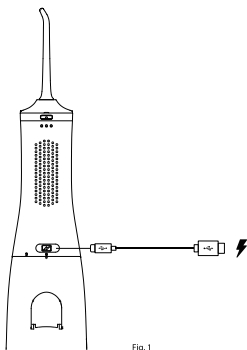
## IRRIGATOR USE

- Place the tip on the body of the device (Fig. 2)
- Open the flap and fill the container with water. The water temperature should not exceed 40°C. Close the container tightly to avoid water leakage (Fig. 3)

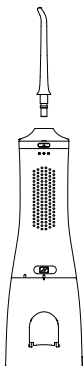
Fig. 2

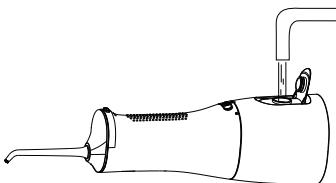
Fig. 3

- Alternatively, you can detach the reservoir from the irrigator and fill it with water through the main opening. To remove the container, turn it counterclockwise (Fig. 4). Once the container is full, attach it back to the body.
- Press the change mode button to select your preferred cleaning programme (Fig. 5). The irrigator stores the last mode used.
  - Soft: gentle cleaning mode, recommended for sensitive gums (2.07-3.45 bar).
  - Standard: standard mode for daily oral hygiene (2.76-5.17 bars).
  - Strengthen: stronger, deeper cleaning of hard-to-reach areas in the mouth (4.48-6.89 bar).
  - Pulse: variable water pressure mode for simultaneous teeth cleaning and gum massage (2.76-7.58 bar).

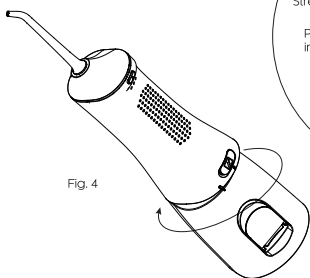
Fig. 4

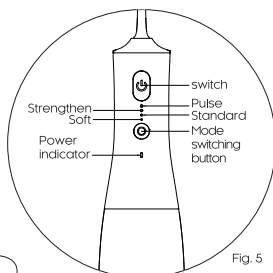
Fig. 5

- Hold the device in an upright position with the nozzle tip pointing towards the teeth.
- Slightly open your mouth to allow the water to flow freely.
- Press the "on/off" button to start irrigation (Fig. 6).
- Direct the water jet at a 90° angle to the teeth and gums.
- Move the irrigator tip slowly along the gum line and interdental spaces.
- Thoroughly clean the areas around braces, crowns and bridges.

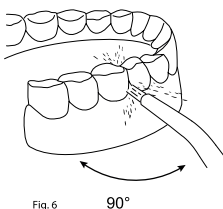
Fig. 6

90°

- You can replace the tip with a different one if necessary. To do this, press the 'on/off' button and then the tip lock release button (Fig. 7).
- When the irrigation is complete, press the 'on/off' button to switch off the device.

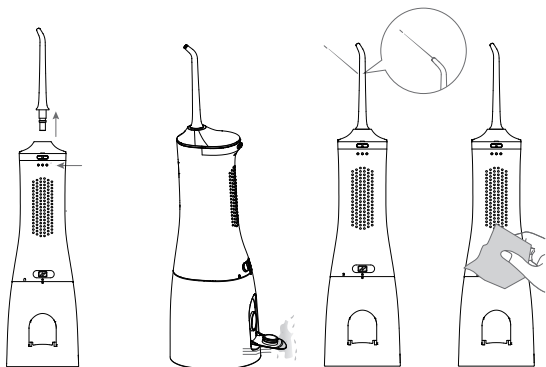
Fig. 7

## AFTER THE IRRIGATION IS COMPLETED (FIG 8.)

- Open the flap of the tank and pour out the remaining water.
- Switch the appliance on briefly using the "on/off" button to rid the appliance of residual water.
- Switch off the unit, close the tank flap and then wipe the units with a dry cloth.
- Use only water and mild cleansers to clean the irrigator. The use of descalers (e.g. for kettles, coffee machines) is not recommended, as they can damage the appliance.

## NOTE

To prevent bacterial growth, check for residual water in the unit. If the appliance is not going to be used for a long period of time, make sure it is properly dried after irrigation.

## TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The irrigator does not work.	The irrigator has just been purchased or has not been used for more than 3 months.	Charge the device for at least 4 hours.
The irrigator only works for a few minutes, even just after it has been charged.	The battery has reached the end of its service life.	Contact the Feelo service centre.
	The charging time was too short.	Charge the device for at least 4 hours.
The water pressure is very low.	The tip is clogged. Hard water promotes limescale build-up, which can clog the tip.	Replace the tip with a new one.
	The tip is deformed.	
Water does not flow out of the tip.	The container is empty.	Fill the container with water.
	The device is tilted too much during use.	Hold the irrigator in an upright position.

## SPECIFICATION

PRODUCT NAME	Feelo PRO Oral hygiene irrigator
Voltage	DC5.0V 1.0A
Power	6.0-7.0 W
Switch type	Button
Display	None
Displayed information	Operation programs + battery charging
Maximum pressure	7.58 bar
Pulsation frequency	1400 pulses/minute
Tank capacity	300 ml
Operation time	Up to 90 minutes on Soft mode
Charging method	Charging via USB-C cable
Charging time	About 4h
Battery type	Li-Ion (2000 mAh)
Noise level	65 dB
Battery-low alert	Flashing charging indicator
Water resistance class	IPX7

## WARNINGS

1. Children and people who are not able to operate the device themselves should not use it. Keep the device out of the reach of children.
2. Never attempt to modify, disassemble or repair the device yourself.
3. Do not place or store the product where it can fall or land in the bath or sink.
4. Do not reach for the product if it has fallen into the water. Immediately disconnect the device from the power supply.
5. Do not use the product while bathing.

## IMPORTANT NOTES

- The irrigator can be used by children aged 8 years and over. Children and persons with limited physical, sensory or mental abilities or lacking experience and knowledge should use the device under the supervision of an adult supervisor, or should be adequately trained in the correct and safe use of the device.
- This product is not a toy. Cleaning and maintenance must not be carried out by children.
- Do not fill the tank with water at a temperature higher than 40°C.
- Do not expose the product to high temperatures or sunlight for a long time. Keep it away from heated surfaces.
- Do not use the device if you lose the nozzle.
- Do not share your tip with other users.
- Do not immerse the device in water.
- The device is designed for cleaning interdental spaces, gums, tongue. Do not use it in any other way.

- Keep the tip away from young children. The nozzle is small enough to be swallowed and stuck in the throat.
- Some types of mouthwash can damage the device, causing the housing or water container to break. Therefore, mouthwash should not be poured into the container.
- If you have undergone oral surgery or have been diagnosed with severe periodontal disease, consult your physician before use.
- Do not press the nozzle too hard against the teeth or gums, as this can damage them.

## TWO-YEAR GUARANTEE

The Feelo water flosser is covered by a 2-year warranty from the date of purchase.

The guarantee covers damage related to manufacturing defects or incorrect workmanship of the product. Keep the proof of purchase for the duration of the guarantee period.

Within the warranty period, the buyer can make a claim to the place of purchase only if the product has been used properly for its intended purpose.

**The guarantee does not cover:**

1. **Surface damage resulting from normal use:**  
Any scratches or abrasions on the surface of the handle shall be regarded as surface damage resulting from normal use and therefore are not covered by the guarantee. The guarantee applies to components that affect irrigator operation.

## 2. **Damage resulting from:**

- Accidental dropping or impact, misuse, abuse or commercial use, negligence or modification of the appliance.
- Improper maintenance, lack of care of the product.
- Connecting to an incorrect power source.
- Repairs or attempted repairs by a party not authorised by Feelo. Any attempts to open or dismantle the appliance, its parts or accessories will render the guarantee invalid. Guarantee claims must be supported by evidence. A claim must be made within the guarantee period.

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