
# Feelo.
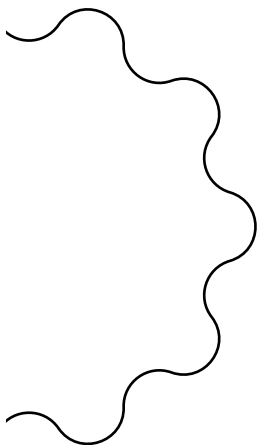
## EASY

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## SET ELEMENTS

- 1 brush handle
- 8 interchangeable tips
- 1 protective cap for the brush head
- 1 travel case
- 1 USB power cable


## BRUSH FEATURES AND FUNCTIONS

- More than 40,000 movements per minute ensure gentle and effective cleaning of the teeth.
- 5 working programs: White (whitening), Clean (daily cleaning), Junior (gentle mode for children over 7 years of age), Polish (polishing), Massage (gum massage).
- Timer – measures 2 minutes, signaling every 30 seconds a break to change the quarter of the dentition.
- Charging via the USB power cable ensures easy and convenient use of the toothbrush at home or on the go.
- Powerful lithium-ion battery – 1 charge is enough for about 25 days of use.
- High degree of IPX7 waterproof. This means that the toothbrush is resistant to accidental submersion in water up to the depth of 1 metre for up to 30 minutes.
- The replaceable brush head has bristles arranged in a wave shape, which ensures thorough cleaning of teeth and hard-to-reach places in the mouth.

## BRUSH DESCRIPTION

1. Brush head
2. On/Off function button
3. Whitening program
4. Program for daily cleaning
5. Program for children from 7 years of age
6. Program for polishing and brightening teeth
7. Gum massage program
8. Charging indicator
9. Charging port

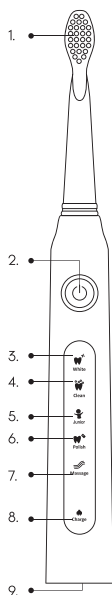
## USING THE BRUSH

- Put the tip on the neck of the handle and then gently twist it clockwise (see: “Fitting and removing the brush head”).
- Wet the head of the tip with water and apply an appropriate amount of toothpaste.
- Place the tip against your teeth at a 45° angle to the gumline.
- Press the On/Off button to turn on the toothbrush. Then keep pressing it to cycle through the 5 modes: White, Clean, Junior, Polish and Massage. The active program is highlighted on the handle panel.
- Gently glide the brush over your teeth and gums. Make sure that the bristles of the brush reach the interdental spaces and massage the gums.
- Do not apply too much pressure to the teeth or apply the brush bristles to one area for a long time.

- The toothbrush turns off after 2 minutes of work, i.e. after the time recommended by dentists.
- To avoid toothpaste splattering, always switch on the toothbrush only after touching it to your teeth, and only remove it from your mouth after brushing your teeth.
- To stop the brush at any time, press the On/Off button.
- With each subsequent use, the toothbrush will turn on on the previously selected program.
- If you want to turn off the toothbrush after selecting the work program, wait 5 seconds. Then press the button to turn off the device.

## FITTING AND REMOVING THE BRUSH HEAD

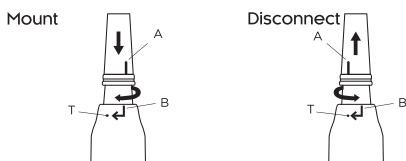- **Mount:** Put the tip on the neck of the brush. Turn it clockwise until point A meets point T.
- **Disconnect:** Turn the tip counterclockwise until point A meets point B. Pull the tip up.

Recommended frequency of replacement of the tip – every 3 months.

## CHARGING THE TOOTHBRUSH

- When the battery is low, the charging indicator flashes 10 times to remind you to charge the toothbrush.
- Remove the plug on the bottom of the handle and connect the toothbrush with the USB cable to your computer or USB charger to start charging.
- When connected, the charging indicator will light up.
- When the battery is fully charged, the indicator will turn off and the charger will turn off automatically.
- When using the toothbrush for the first time or if the toothbrush has not been used for a long time, it may take up to 5 minutes for the charging indicator to light up.
- Charge your toothbrush only when it needs to. Continuous charging of the device is not recommended.

**ATTENTION:** Do not use a damaged cable or charger for charging.

## CLEANING AND MAINTENANCE

- After brushing your teeth, rinse the toothbrush under running water.
- Remove the tip and also rinse it under running water. Also rinse the neck of the handle. Shake the tip to get rid of any residual water. Once a week, rinse the entire tip, including the bottom, with warm water.
- Use a clean, slightly damp cloth to clean the entire surface of the handle. Wipe the handle dry.

## SOLVING PROBLEMS

PROBLEM	SOLUTION
The toothbrush won't turn on.	Charge your toothbrush.
	If the toothbrush still does not start after the charging, contact the customer service department.
The vibration of the brush decreases.	Charge your toothbrush or replace the brush head. The brush vibrates differently depending on the selected work program.
The charging light does not light up when the toothbrush is connected to charge.	If the toothbrush has not been used for a long time, it enters a sleep mode. After 5 minutes, the toothbrush should awaken from the sleep mode and begin charging the battery. This is a standard phenomena.
	If the toothbrush still does not charge, contact the customer service department.
Tickling when brushing teeth.	The movements of the sonic toothbrush can cause a tickling sensation, especially the first time you use it. This impression will disappear with time.

## SPECIFICATION

PRODUCT NAME	Feelo Easy Sonic Toothbrush
Rated power	1 W
Input	DC5V
Switch type	Button
Display	None
Information displayed	Work programs + battery charging
Operation time	Up to 25 days
Way of charging	Charging via USB cable
Charging time	About 4h
Battery	260 mAh
Battery discharge signalling	The charging indicator flashes
Water resistance class	IPX7

## WARNINGS

1. Do not disassemble or attempt to repair the product yourself.
2. If you have undergone oral surgery or have been diagnosed with severe periodontal disease, consult your physician before use.
3. Children and persons unable to operate the product on their own should only use it under appropriate adult supervision.
4. Do not immerse the product in water.
5. Do not expose the product to high temperatures or sunlight for a long time.
6. Do not press the toothbrush hard against the gums and teeth.

## TWO-YEAR GUARANTEE

The Feelo toothbrush is covered by a 2-year warranty from the date of purchase. The guarantee covers damage related to manufacturing defects or incorrect workmanship of the product. Keep the proof of purchase for the duration of the guarantee period. Within the warranty period, the buyer can make a claim to the place of purchase only if the product has been used properly for its intended purpose.

**The guarantee does not cover:**

1. **Surface damage resulting from normal use:** Scratches or abrasions on the surface of the of the handle will be treated as surface damage resulting from normal use and will therefore not be covered by the warranty. The guarantee applies to the components that affect the functioning of the toothbrush.

## 2. **Damage resulting from:**

- Accidental dropping or impact, misuse, abuse or commercial use, negligence or modification of the appliance.
- Improper maintenance, lack of care of the product.
- Connecting to an incorrect power source.
- Repairs or attempted repairs by a party not authorised by Feelo. Any attempts to open or dismantle the appliance, its parts or accessories will render the guarantee invalid. Guarantee claims must be supported by evidence. A claim must be made within the guarantee period.

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