
Feelo.

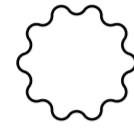KIDS

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EN

SET ELEMENTS

- 1 brush handle
- 3 replaceable heads 3-6 years old
- 1 AAA battery
- Instruction manual

BRUSH FEATURES AND FUNCTIONS

- The toothbrush makes 24,000 vibrations per minute, providing gentle and effective cleaning of your child's teeth and gums.
- Colourful LED illumination of the toothbrush head encourages children to brush their teeth.
- Ergonomic handle with non-slip coating made from ABS material.
- Timer - 2-minute automatic timer, reminding every 30 seconds to move on to the next quadrant of the mouth. The toothbrush turns off automatically after two minutes.
- Power Source: 1 AAA battery. The average operating time of a battery is approximately 150 days when using the toothbrush 2 times a day for up to 2 minutes each time.
- IPX7 water-resistant rating. This means that the toothbrush is resistant to accidental submersion in water up to the depth of 1 metre for up to 30 minutes.
- The replaceable toothbrush head is gentle to gums thanks to its DuPont bristles.

BRUSH DESCRIPTION

1. Brush head
2. LED head illumination
3. On/Off Button
4. Battery cover

USING THE BRUSH

- Put the head on the holding pin of the handle and then gently twist it clockwise (see: 'Fitting and removing the brush head').
- Wet the head of the head with water and apply an appropriate amount of toothpaste.
- Place the head against your teeth at a 45° angle to the gum line.
- Press the On/Off button. Toothbrush head illuminates white. White LED keeps lighting for 30 seconds.
- Then push the button for a second time to turn on the brush strokes. The lights on the head turns on. You can change the light colour (white or colourful) by pressing the button.
- Start brushing your teeth. Use gentle strokes to brush your teeth and gums; use gentle sweeping motions. Make sure that the bristles of the brush reach the interdental spaces and massage the gums.

EN

- Remember to thoroughly clean the outside, inside and chewing surfaces.
- Timer uses pulse vibrations every 30 seconds to remind to move on to the next quadrant of the mouth.
- Do not apply too much pressure to the teeth or apply the brush bristles to one area for a long time.
- The toothbrush turns off after 2 minutes of work, i.e. after the time recommended by dentists.
- After toothbrushing rinse the toothbrush with water and leave it to dry.
- To avoid toothpaste splattering, always switch on the toothbrush only after you have placed it against your teeth, and only remove it from your mouth after brushing your teeth.
- To stop the brush at any time, press the On/Off button.

FITTING AND REMOVING THE BRUSH HEAD


- **Fitting:** Put the head on the holding pin of the handle. Gently twist it clockwise. Make sure that the notch on the holding pin fits into the corresponding groove inside the replaceable head.
- **Removing:** Twist the head counterclockwise. Pull the head upwards.

Recommended frequency of replacement of the head - every 3 months.

REPLACING THE BATTERY

- Looking from the back of the toothbrush, twist the battery cover clockwise ①, so that the line indicator is aligned with the left end of the arrow. Then pull the cover downwards ②.
- Insert a new AAA battery remembering about the correct polarity.
- Put the lid on ③, so that the line indicator is aligned with the left end of the arrow. Then twist it counterclockwise so that the line indicator is aligned with the dot ④.

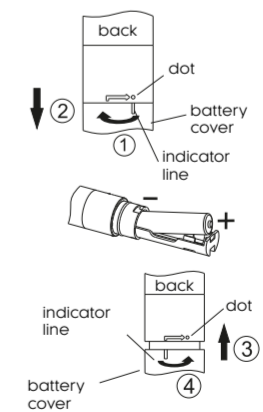EN

CLEANING AND MAINTENANCE

- After toothbrushing remove the head from the handle and rinse both parts with water.
- Shake the head to get rid of any residual water.
- Use a clean, slightly damp cloth to clean the entire surface of the handle. Wipe the handle dry.

SOLVING PROBLEMS

PROBLEM	SOLUTION
The toothbrush doesn't turn on.	Replace the battery. If the toothbrush still does not turn on after replacing the battery, contact the customer service.
The strokes are getting weaker.	Replace the battery or the head. The strokes may get weaker if there is too much pressure applied.
Short battery life.	Replace the battery. If the toothbrush still has short battery life after replacing the battery, contact the customer service.

SPECIFICATION

PRODUCT NAME	Feelo Kids Sonic Toothbrush
Material	DuPont Nylon/ABS/TPE
Power Source	1 AAA battery
One battery life	Up to 150 days
Brushing strokes frequency	24,000 strokes per minute
Water resistance class	IPX7
Dimensions	188 x 19 x 19 mm
Quality standard	GB4706.59-2008

WARNINGS

- This product is not a toy. Children and persons unable to operate the product on their own should only use it under appropriate adult supervision.
- Do not disassemble or attempt to repair the product on your own. If your child has recently undergone a surgery within the oral cavity, please consult your doctor before using this product.
- Remove the battery if you do not plan on using the device for a longer period of time.
- Insert the battery whilst paying attention to the correct polarity. Charge only batteries that are designed to be charged. Rechargeable battery should be removed before charging.
- The toothbrush should be used as intended, namely for toothbrushing and cleaning of gums and tongue. The toothbrush shall not be used for any other purposes.
- Do not immerse the product in water.
- Do not expose the product to high temperatures or sunlight for a long time.
- Do not press the toothbrush hard against the gums and teeth.

TWO-YEAR GUARANTEE

The Feelo toothbrush is covered by a 2-year warranty from the date of purchase. The guarantee covers damage related to manufacturing defects or incorrect workmanship of the product. Keep the proof of purchase for the duration of the guarantee period. Within the warranty period, the buyer can make a claim to the place of purchase only if the product has been used properly for its intended purpose.

The guarantee does not cover:

- Surface damage resulting from normal use:** Any scratches or abrasions on the surface of the handle shall be regarded as surface damage resulting from normal use and therefore are not covered by the guarantee. The guarantee applies to components that affect the functioning of the toothbrush.
- Damage resulting from:**
 - Accidental dropping or impact, misuse, abuse or commercial use, negligence or modification of the appliance.
 - Improper maintenance, lack of care of the product.
 - Repairs or attempted repairs by a party not authorised by Feelo. Any attempts to open or dismantle the appliance, its parts or accessories shall render the guarantee invalid. Guarantee claims must be supported by evidence. A claim must be made within the guarantee period.

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