# S·E·Y·S·S·O

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- 2 brush heads
- 1 handle
- 1 USB cable

### PLUGGING AND CHARGING

- The toothbrush is charged using a USB-DC 5V cable.
- The first full charge of the toothbrush takes about 4 hours and lasts up to 3 weeks (on the Soft program), which is about 42 brushing cycles.
- Each subsequent charge of the toothbrush should take about 4 hours.
- During charging, you must disconnect the toothbrush from the cable for the time of brushing your teeth.
- After a two-minute tooth brushing time, you can reconnect the toothbrush to charge again.
- There is no risk of overcharging the battery.
- To extend the battery life, it is recommended to charge the toothbrush when it signals a low battery level (flashing of three operation LEDs in white).
- After the battery is fully charged, all LEDs light up briefly and then go off.

## USING THE SEYSSO DINOSAUR SONIC BRUSH

Place the head on the handle. Then apply the appropriate amount of toothpaste to the bristles. Before turning on the toothbrush, place the head in the mouth. Turn on the toothbrush. After the appropriate time, depending on the program, the device will turn off automatically and is ready for reuse. During brushing, through an interval that occurs every 30 seconds, the toothbrush informs us of the need to change the tooth brushing zone. The toothbrush moves at a speed of 70,000 sonic movements per minute in three intensity settings, within the Clean, Soft, and Massage programs.

## HOW TO CHANGE THE BRUSHING MODES?

The toothbrush's power button also controls the selection of the operating program. It allows us to change the operating speed sequentially between: Clean, Soft, and Massage. Use the button to turn on the device, then press it further to switch between the toothbrush's programs. If a program is active for more than 5 seconds, the next press of the button will turn the device off. The SEYSSO Dinosaur sonic toothbrush turns on with the Clean program first.

- To switch from the Clean program to the next, press the function selection button a second time.
- To change to the next operating program, continue pressing the function button.
- After 5 seconds, the working program is memorized and the next use of the button will turn the toothbrush off.
- When the toothbrush is turned on again, it will operate on the last remembered program.

## WARNINGS! READ BEFORE USING THE TOOTHBRUSH!

SAVE THIS MANUAL FOR FURTHER REFERENCE.

- For home use only. 
- The toothbrush is equipped with a waterproof handle. It is safe and designed for bathroom use. The device has an IPX7 protection rating. This means that the toothbrush is resistant to accidental immersion in water up to 1 meter deep for up to 30 minutes.
- Remember: the toothbrush is not adapted to operate under water!
- The product is not intended for self-repair.
- Straighten the cable before use.
- In case of a failure, do not disassemble the device and do not replace any components, as there is a risk of electric shock. Contact the manufacturer's service.
- A damaged power cord must be replaced with a new one. Contact the manufacturer's service.
- If the toothbrush is dropped, replace the head before the next use. Even if the damage is not visible.

## MAINTENANCE AND CLEANING

- After brushing your teeth, rinse the toothbrush under running water.
- Remove the head from the handle and rinse it under running water.
- Do not use water hotter than 40°C or acidic cleaning solutions.
- Dry the handle before connecting it to the charger.
- Before cleaning the toothbrush, make sure they are disconnected from the power source. Use only a dry cloth to wipe the surface of the charger! With any electrical device, especially in the presence of children, observe basic precautions.

## BRUSH HEAD FITTING AND REMOVAL

Before using the toothbrush for the first time, apply a thin layer of lubricant onto the brush holder metal pin. This will make it easier to fit and remove the toothbrush heads. We recommend changing brush these heads every 3 months.

### WARNING - DANGER

To reduce the risk of electric shock:

- Always disconnect the device from the power source after use.
- Do not use while bathing.
- For more detailed instructions on storage and use of the device, as well as included accessories, refer to SEYSSO.

To limit the risk of fire, short circuit, or injury:

- Disconnect the charger from the power supply after charging is complete.
- The product is not a toy. Exercise particular caution if the toothbrush is used by a child or in their presence.
- Use as intended. Use only accessories recommended by the manufacturer.
- Do not use the charging cord if it has been damaged. If after dropping, leaving in the open air, or submerging in water the toothbrush does not work properly, return it to the place of purchase.
- Do not stretch the charging cord, cut it with doors, wrap it around sharp edges. Do not turn on the toothbrush while charging. Keep the components of the device away from heated surfaces.

- Do not hold the plug with wet hands.
- Never use the product if you have problems with concentration and feel drowsy.
- Do not use the device outdoors or in places where it may be exposed to aerosols or pure oxygen.
- Before connecting the device to the socket, familiarize yourself with its specification.

It is not recommended to use the device by children under the age of 3 years. Children and people with limited sensory, motor, or mental abilities should use the device under the supervision of a guardian or should be properly trained in the correct and safe use of the device. Children without supervision should not perform cleaning and maintenance of the device.

Do not use too much force or too fast movements while brushing your teeth. In the first few days of using the toothbrush, the gums may bleed slightly. The bleeding should stop on its own after a few days. If it persists or increases after this time, it is recommended to contact a dentist. For delicate gums and sensitive teeth, it is recommended to use the Soft program.

## ENVIRONMENTAL PROTECTION

Do not dispose of the SEYSSO Dinosaur toothbrush with municipal waste. Take it to an electronic equipment disposal and recycling point.

## TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
Toothbrush movement seems weaker.	Pressing the toothbrush too hard against the teeth.	Reduce the pressure of the toothbrush.
After loading the toothbrush, the battery life is very short.	Incorrect positioning of the toothbrush on the charger can make the charging process ineffective.	Make sure that the toothbrush is placed correctly when charging.
		If the problem persists, battery wear is probably the cause. Contact the customer service department.
		If the toothbrush still does not charge, contact the customer service department.
The toothbrush does not start.	Check if it is caused by insufficient supply power.	Check if the toothbrush has been charged correctly.
		If the toothbrush still does not start after charging, contact the customer service department.

## SEYSSO DINOSAUR TOOTHBRUSH SPECIFICATIONS

PRODUCT NAME	SEYSSO DINOSAUR
Switch type	Button
Display	LED diode
Operation time	Up to 84 minutes
Charging method	USB-DC 5V cable charging
Charging me	Approx. 4 h
Input rated voltage	5 V
Input rated power	1,5 W
Battery capacity	600 mAh
Battery	li-ion
Battery-low alert	LED diode
Water resistance class	IPX7

Important: Use of the correct power adapter. Please ensure that the power adapter used has output parameters that are in accordance with the parameters listed on the label of your toothbrush. Using a power adapter with incompatible parameters can lead to improper functioning of the device and is absolutely unacceptable. Failure to adhere to this rule may result in damage to the device and may also lead to the loss of warranty. Take care of your device by using only recommended accessories and power adapters.

## TWO-YEAR GUARANTEE

SEYSSO provides a guarantee for the period of 2 years from the date of purchase. The guarantee covers damage related to manufacturing defects or incorrect workmanship of the product. Keep the proof of purchase for the duration of the guarantee period. A complaint can be filed during the guarantee period if the product has been used correctly and in line with its intended use.

### THE GUARANTEE DOES NOT COVER:

1. Surface damage resulting from normal use.
2. Damage resulting from:
  - a) accidental dropping or impact, misuse, abuse or commercial use, negligence or modification of the appliance,
  - b) improper maintenance, lack of care of the product,
  - c) connecting to an incorrect power source
  - d) repairs or attempted repairs by a party not authorised by SEYSSO

Any attempts to open or dismantle the appliance, its parts or accessories will render the guarantee invalid. Guarantee claims must be supported by evidence. A claim must be made within the guarantee period. The material covering the brush holder is delicate, so it is recommended to exercise caution during use. Scratches or abrasions on the surface of the brush holder will be treated as surface damage resulting from normal use. Thus, they will not be covered by the guarantee. The guarantee applies to components that affect toothbrush operation.

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Water Resistance (Toothbrush Handle): IPX7. Protection against the effects of immersion in water up to a depth of 1 meter for 30 minutes.  
The toothbrushes are equipped with a built-in lithium-ion battery.

 This product meets the requirements of the so-called "New Approach" Directives of the European Union (EU).

The packaging is recyclable.

The instruction manual is recyclable.


This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU with the symbol of a crossed-out waste bin. This marking indicates that the equipment, after its period of use, should not be disposed of with other household waste. The user is obliged to hand it over to entities collecting waste electrical and electronic equipment, including local collection points, stores, and municipal units, which create an appropriate system allowing for the disposal of this equipment. Proper handling of waste electrical and electronic equipment helps avoid harmful consequences for human health and the environment due to the presence of hazardous components and improper storage and processing of such equipment. In the device, the occurrence of certain substances that can negatively affect the environment during the use of this product and after its disposal has been limited.

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Creative Notebook  
with coloring pages  
and puzzles

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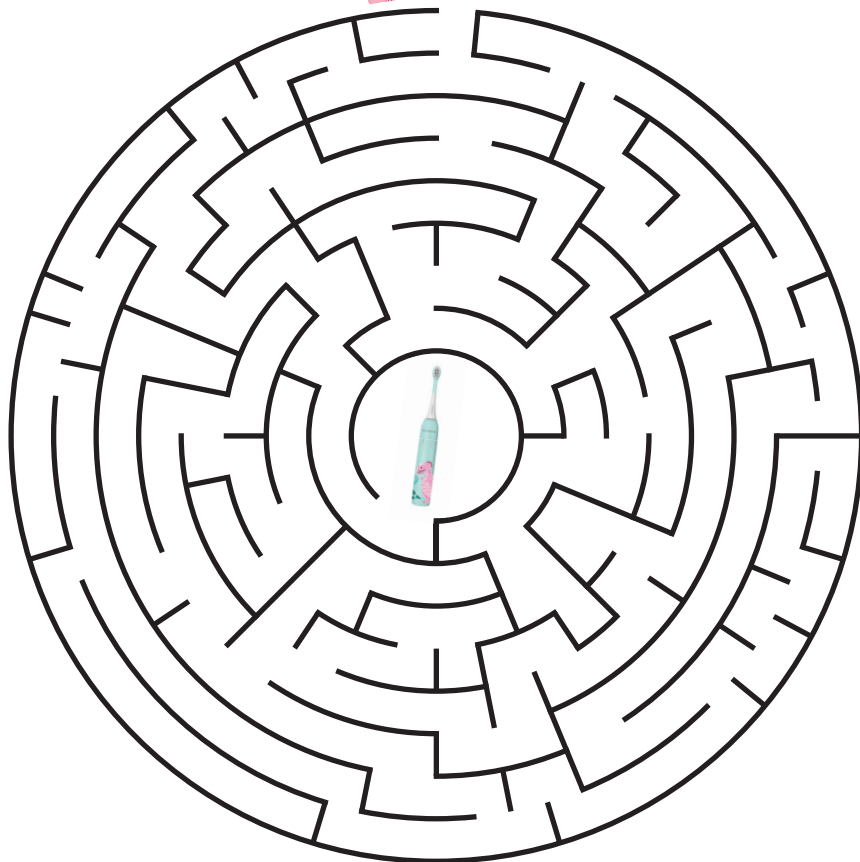
## Learn 10 rules for a healthy smile:

1. I brush my teeth for at least 2 minutes.
  2. I do it at least twice a day.
  3. I apply the appropriate amount of toothpaste to the brush.
  4. I use a soft toothbrush to brush my teeth.
  5. I don't swallow the toothpaste.
  6. After brushing my teeth, I spit out the toothpaste and do not rinse my mouth with water.
  7. I don't chew candies and lollipops.
  8. I don't drink sweet drinks in the evening.
  9. I don't brush my teeth immediately after eating sour things.
  10. I visit the dentist at least twice a year.
- 

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
Help seyssosaurus find his lost toothbrush.


pokoloruj obrazek.

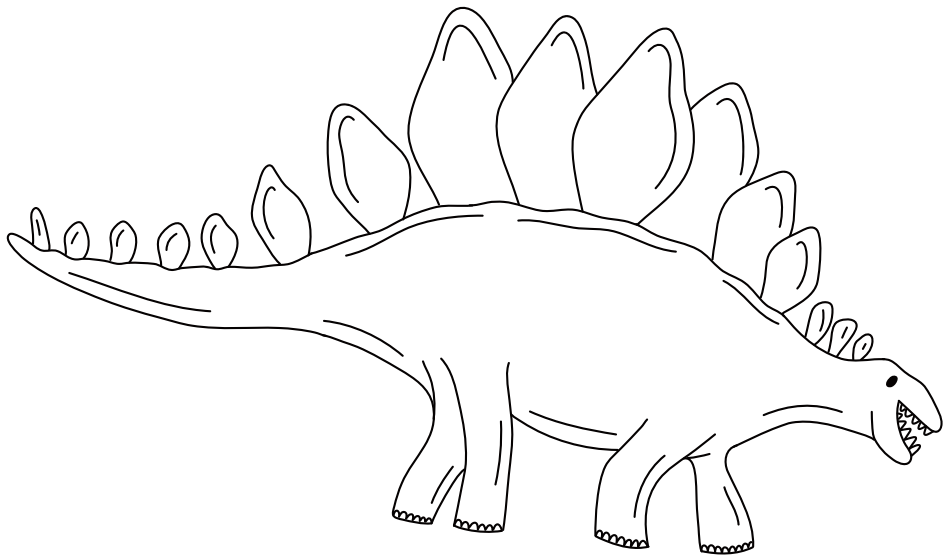
color the picture.

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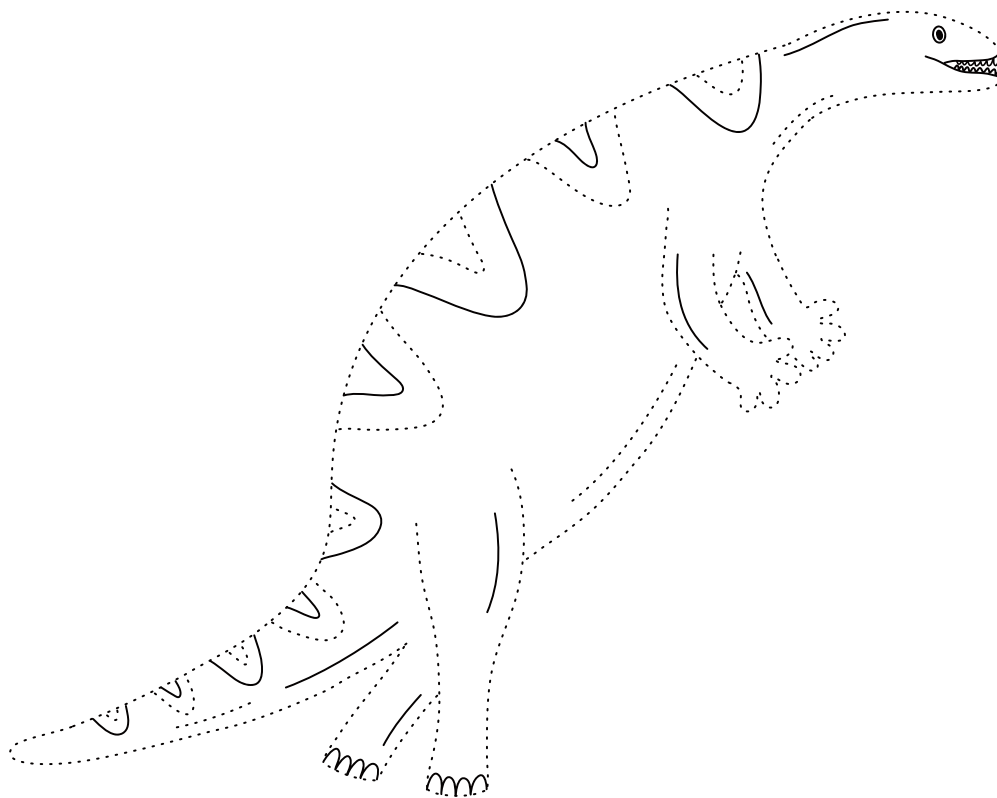
connect the dots and color the picture.

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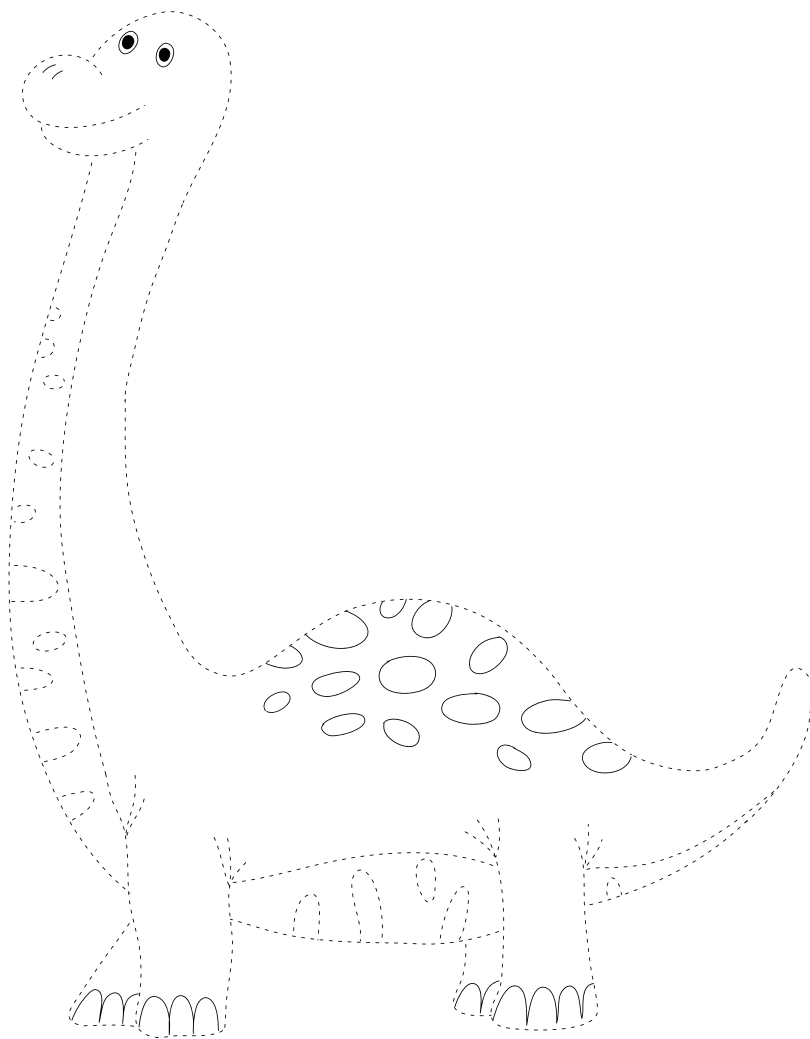connect the dots and color the picture.

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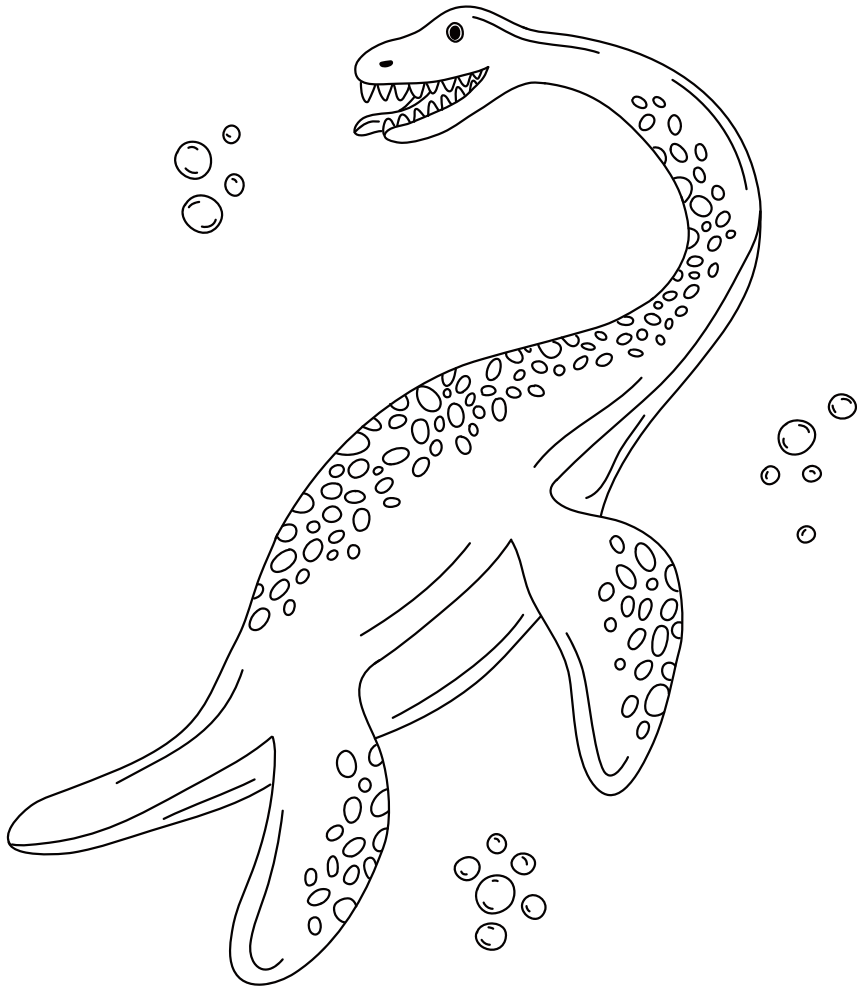
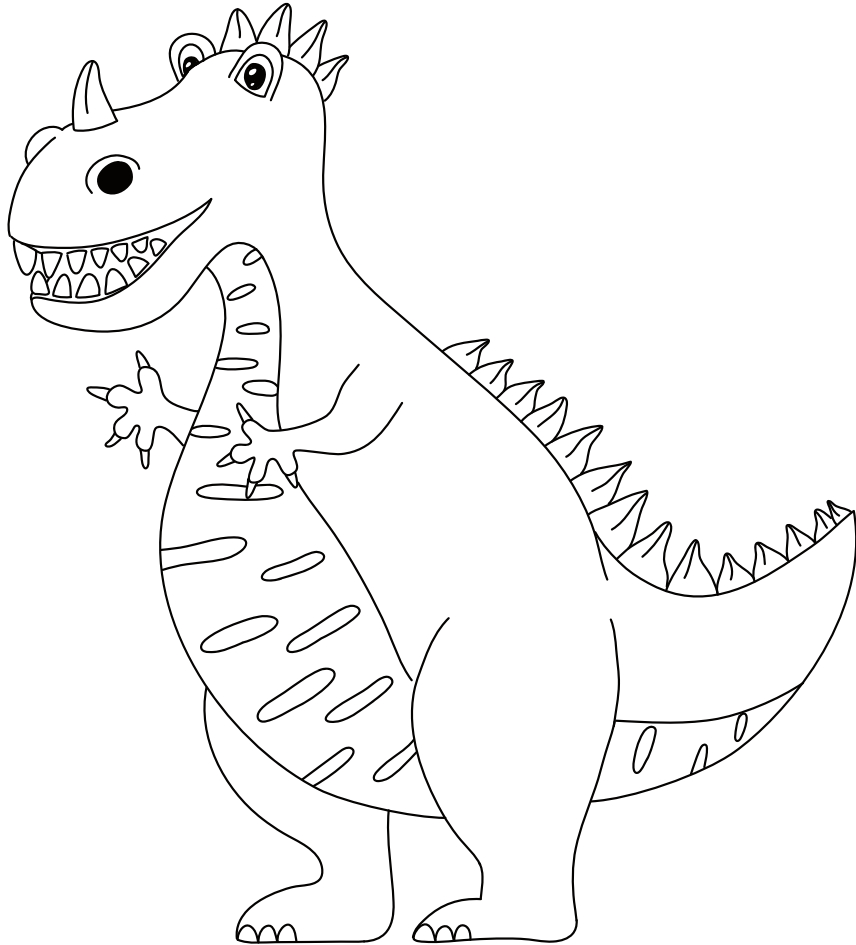
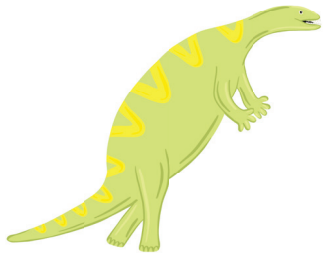
match the dinosaurs with their toothbrushes.

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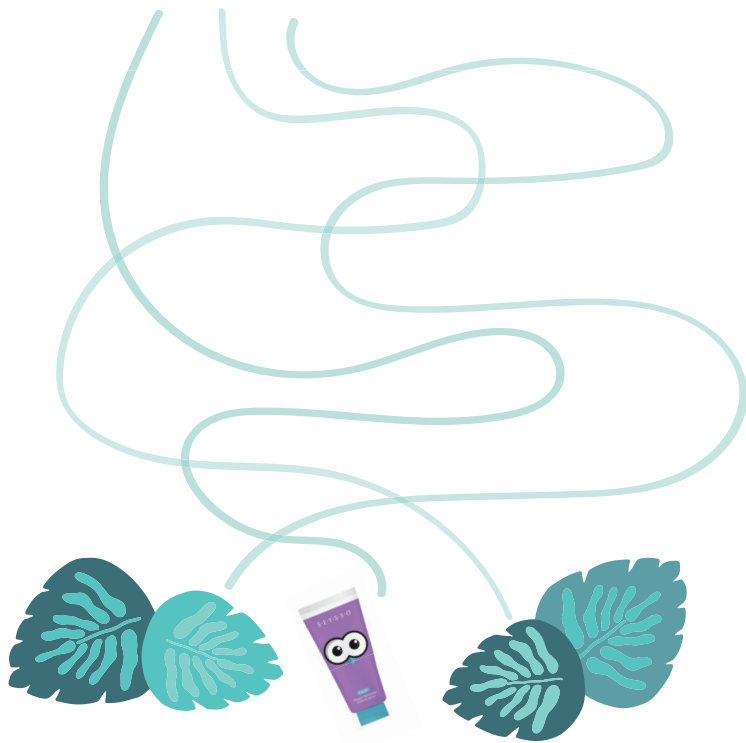
Which path will lead the dinosaur to the toothpaste?

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



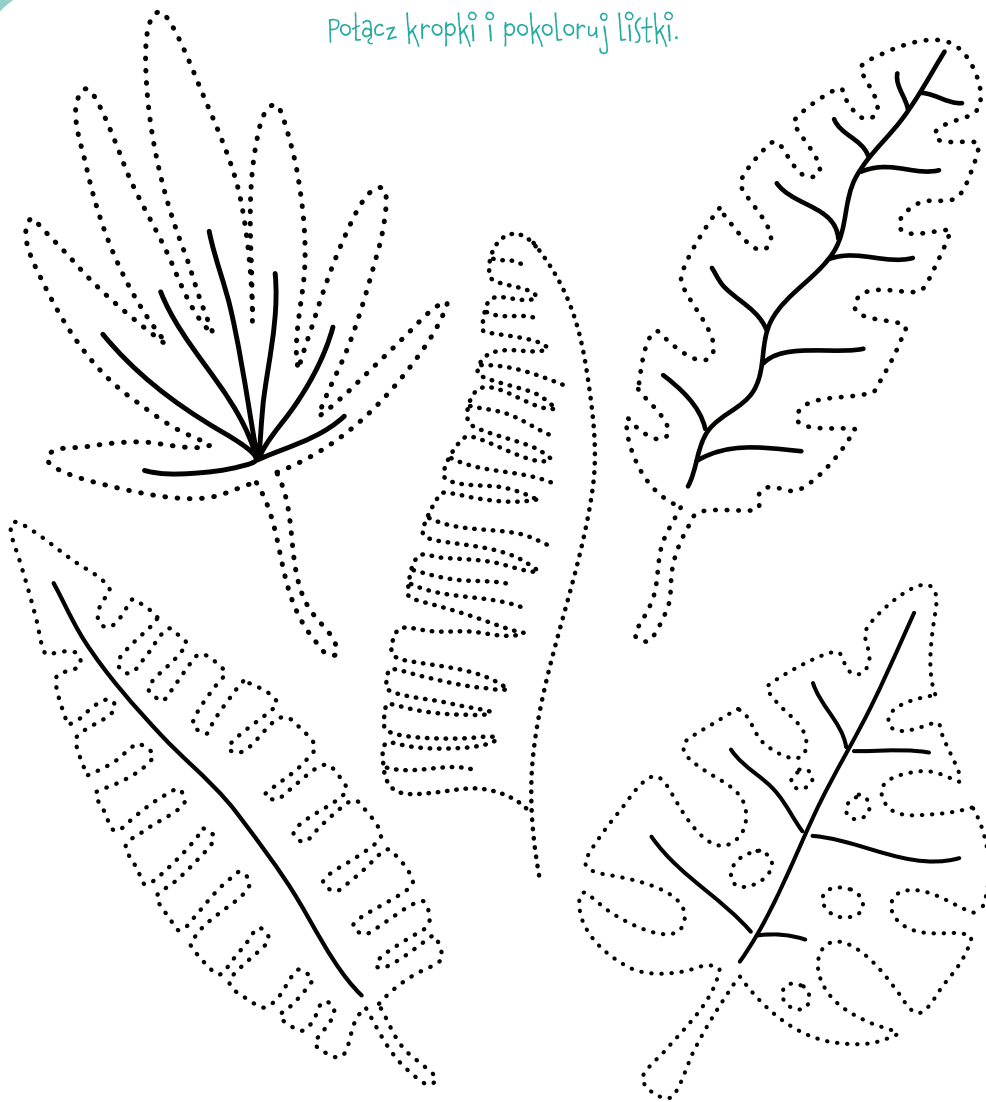
connect the dots and color the leaves.

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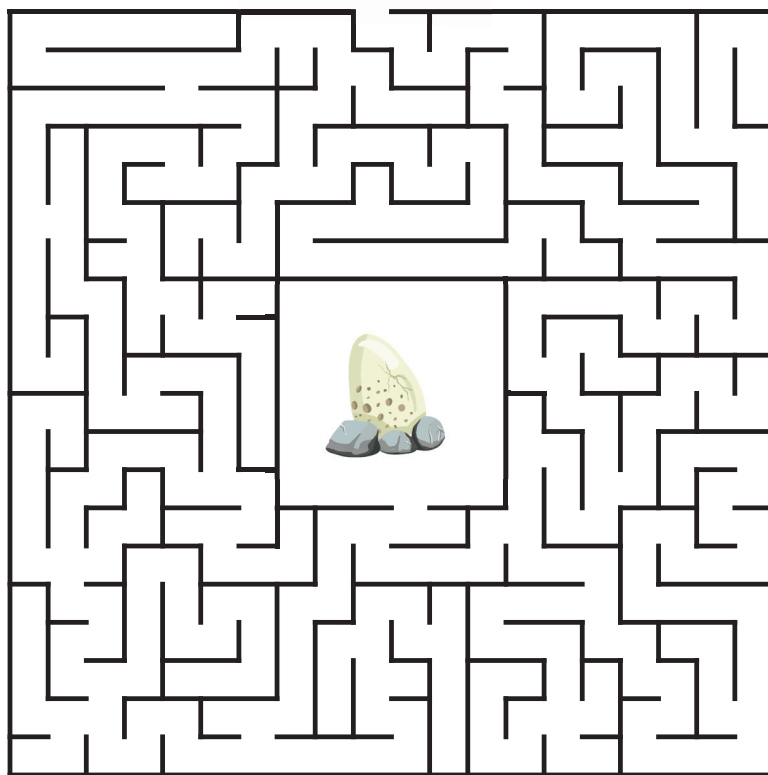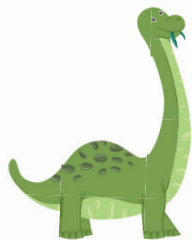Help the dinosaur find the way to the egg.

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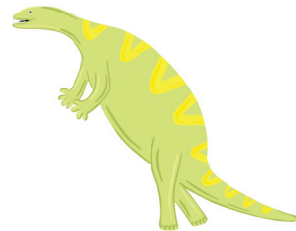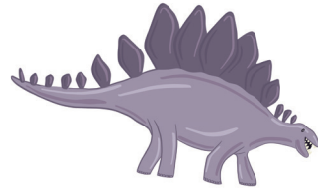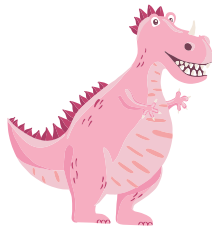Do you remember what these dinosaurs are called?