S · E · Y · S · S · O

GOLD FOREST GREEN
SONIC TOOTHBRUSH

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
User Manual

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Before using this product, please read carefully this user's manual and keep it properly for later reference.

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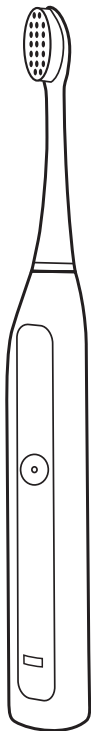English

1	Attention	4
2	Articles in the packaging box	7
3	How to charge your SEYSSO sonic toothbrush	8
4	Disinfection of the SEYSSO heads with UV light in the travel case	9
5	Charging the SEYSSO brush in the travel case	10
6	Washing the teeth with the SEYSSO brush step by step	10
7	Washing the tongue with the SEYSSO brush step by step	11
8	Cleaning and massaging the skin of the face with the SEYSSO brush step by step	12
9	How to clean your SEYSSO sonic toothbrush	13
10	How to store your SEYSSO sonic toothbrush	14
11	Environmental protection	14
12	Replacement heads	14
13	Warranty and services	15

1 Attention

Before using this product, please read carefully this user's manual and keep it properly for later reference

Danger:

- Keep the inserted charging base away from water. Do not place it above or nearby any container with water (e.g. bathtubs or washbasins). Do not place the charging base or rechargeable travel case in water or any other liquid. After cleaning the product and before connecting it to the power, please ensure that the charging base or the rechargeable travel case is thoroughly dry. Do not use the charging station or the case if they have fallen into the water.

Warning:

- Before you connect the charging station or the travel case to the mains, make sure that the voltage of the electrical socket corresponds to the requirements of the appliance.
- If the charging base or the travel case has any damage, please stop using the product. If the charging base, travel case, power line or plug is damaged, please be sure to use a new one of the original model to replace the damaged one to avoid danger.
- If any part of the kit, such as the replacement head or the brush handle, is damaged, do not use it. Do not attempt to repair, modify, or disassemble the kit parts yourself. In the event of damage, contact the customer service center in your country or region (see the „Warranty and Services“ section).
- Do not plug the travel case to the mains if the UV light bulb does not switch off when the case is open. UV light can be harmful to the eyes and skin.
- Do not plug the travel case to the mains if it is impossible to close it properly due to damage.
- Do not connect the case to the mains in places where aerosols are sprayed or oxygen is administered.

- Please do not use charging base or travel case in the open air or nearby an article with high-temperature surface.
- Do not charge, use or leave the brush exposed to high temperatures or external influences.
- The brush, charging station and case are not intended for use by unattended children or persons with reduced motor, mental or cognitive abilities. Keep all components of the kit out of their reach. The product is not a toy.

Attention:

- If you had oral or gum operation in the past two months, please consult your doctor before using this product.
- The brush may cause slight bleeding of the gums during the first days of use. If bleeding continues within a week or is extremely heavy, consult a physician.
- The product complies with all standards for the use of electromagnetic appliances. If you have any questions about the use of the brush by persons with cardiac devices.
- The brush is intended to clean the teeth, gums, tongue and the skin of the face (depending on the head selected). Do not use it in any other way. If you experience pain or discomfort after using the brush or have any concerns about the health effects of the brush, consult a physician.
- If you have any medical problems, please consult your doctor before using this SEYSSO sonic toothbrush.
- Please stop using deformed brush heads or those with curved bristles. Please replace your brush head every three months. If the brush head is worn, please replace it ahead of time.
- Please do not clean your brush head in a dish-washing machine.
- Use the brush, travel case, charging station and other components only for their intended purpose.
- Please do not use your brush heads not recommended by the manufacturer.

- If the toothpaste you use contains peroxide, sodium bicarbonate or other bicarbonate (commonly seen ingredients in whitening toothpaste), it is suggested that you clean your brush head thoroughly with soap-suds after using it every time.


EMF:

This product accords with all relevant standards for EMF. According to current scientific proof, if you use this product correctly in accordance with this user's manual, this product is safe.

UV-C light:

The travel case is equipped with a UV-C light bulb emitting waves of 253.7 nm in length and 700-1500 $\mu\text{W}/\text{cm}^2$ in intensity. Direct contact with UV-C light can have adverse effects on your health. The tightly closed case fully protects against UV-C light. The bulb is non-replaceable, do not remove it.

2 Articles in the packaging box

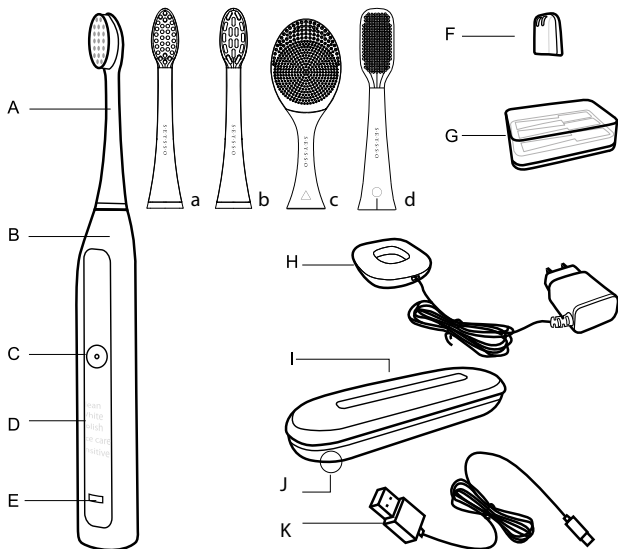- A. Brush head
- a) SEYSSO Gold Black Anti-plaque 2 pcs.
- b) SEYSSO Gold Black Whitening 1 pc.
- c) SEYSSO Gold Black Dermo Sonic 1 pc.
- d) SEYSSO Gold Black Tongue Cleaner 1 pc.
- B. Brush handle
- C. Power switch button
- D. Tooth brushing modes
- E. Charging indicator
- F. Travel brush head protective cap 3 pcs.

- G. Brush heads storage box
- H. Inductive charging base
- I. Travel case with the function of brush charging and brush head disinfection
- J. Case LED
- K. MicroUSB charger


3 How to charge your SEYSSO sonic toothbrush

Insert the toothbrush handle into the charging base to charge it. Electricity quality is shown on the LED electricity quality display. The fast blinking red indicator means low battery (for less than three times). The charging indicator should blink through the charging process and turn to continuous light when battery is fully charged. It will take about 15 hours to fully charge the handle, and it can be used up to 70 days straight.

Result depends on a chosen brushing mode and regards using twice a day:

- a. **Clean**:>57 days;
- b. **White**:> 25 days;
- c. **Polish**:> 45 days;
- d. **Face Care**:> 58 days;
- e. **Sensitive**:> 70 days;

If the battery is low, the electricity quality display LED light will flash and the speed of motor will reduce. When the battery runs out, the motor will stop working. It needs at least 30 minutes to change for reboot and hour later it can start to work. For convenient use, it is suggested that you use this product after it is fully charged.

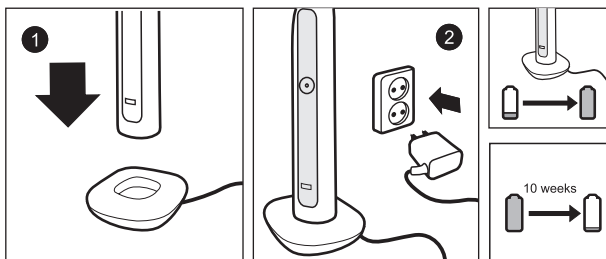4 Disinfection of the SEYSSO heads with UV light in the travel case

The SEYSSO Gold Forest Green travel case has the toothbrush head disinfection function. The case is only compatible with SEYSSO Gold Anti-plaque and SEYSSO Gold Whitening heads. UV-C light emitted by the bulb removes bacteria, viruses, and fungi on the bristles of the toothbrush head.

Note: Direct exposure to UV light is harmful to the eyes and skin.

- Only one head can be disinfected at a time. The UV light bulb is located near the tip of the head on the left side of the case.
- Insert the head into the case. Make sure that it is inserted correctly, without the protective cover, with bristles toward the UV light bulb

Note: if you want to disinfect the head immediately after use, rinse it under the tap and shake off excess water from the head bristles.

- Close the case and plug it into the electrical socket via the microUSB port.

Note! Do not use the case if it is impossible to close it properly due to damage.

- Full disinfection of the head takes 10 minutes. The case LED emits continuous blue light during disinfection. Opening the case will terminate the disinfection. Disinfection of less than 10 minutes may be ineffective.

Note! Do not use the case if the UV light bulb does not switch off when the case is open. UV light can be harmful to the eyes and skin.

- If the brush handle is inside the case, it will charge at the same time as the head is disinfected. When the disinfection is complete, the LED will change from blue to red and the case will automatically switch to brush charging mode.
- If there is no brush inside the case, the case LED will turn off after disinfection is completed.


Note! UV disinfection of the heads may cause the color of the bristles to change gradually.

5 Charging the SEYSSO brush in the travel case

The SEYSSO travel case has the function of brush charging. Put the brush in the case and plug it into the electrical socket via the microUSB port.

Close the case. Opening the case will discontinue the charging. During the first 10 minutes, the brush charging will take place simultaneously with the disinfection of the toothbrush head (see item 4). During this time, the case LED emits continuous blue light. Disinfection is performed whether or not there is a head inside the case.

After 10 minutes, the disinfection of the heads ends and the LED turns red. A blinking red LED indicates that the brush is charging. If the LED emits continuous red light, the brush has been fully charged.

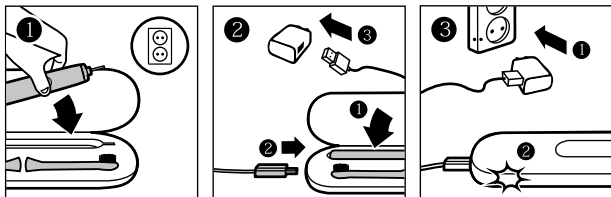
Note! The charger adapter is not included.

6 Washing the teeth with the SEYSSO brush step by step

- Slide the head onto the handle. Then apply the appropriate amount of toothpaste to the brush bristles. We recommend using a sonic-friendly toothpaste. We recommend the SEYSSO toothpaste, which makes full use of the potential of the sonic brushes thanks to its ideal composition and consistency.
- Place the brush head in the mouth before you switch on the brush. Position the toothbrush so that its bristles are at a 45-degree angle toward the

gum line and switch on the toothbrush by pressing the function button. To change the program, press the function key again.

- Move the brush slowly along the gum line but do not press it firmly against the teeth. The toothbrush has a built-in automatic timer that helps to comply with the brushing time recommended by dentists. With a 30-second interval, indicated by a pause in the toothbrush operation, you can evenly divide the brushing time into all four teeth quarters.
- The smart memory chip will remember your preferences and the next time you turn on the toothbrush, it will set the program you selected previously.

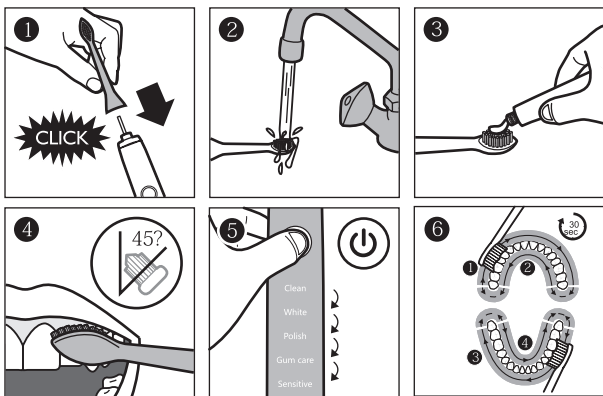7 Washing the tongue with the SEYSSO brush step by step


- Attach the tongue brush to the SEYSSO Gold brush handle.

- Place the brush head on your tongue and turn your toothbrush on with the Clean program selected.
- Massage the tongue, the inner part of cheeks and palate with the silicone nubs of the brush.
- After using the brush, rinse your mouth and brush with lukewarm water.

8

Cleaning and massaging the skin of the face with the SEYSSO brush step by step

- Attach the face brush to the SEYSSO Gold brush handle.
- Put a small amount of gel or face cleansing foam on silicone nubs or directly on your face.
- Wet your face with water and choose the Face Care program.
- Move the brush in a circular motion, changing the side of the tip. Cleansing the eye area, cheeks and neck choose side of the tip with finer silicone nubs. For the T zone, choose side with thicker ones.
- When brushing is complete, rinse your face and the brush with lukewarm water.
- Apply a small amount of cream on a dry face and pat it gently with your fingers.

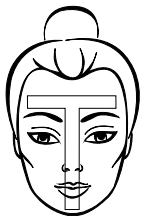Five modes of SEYSSO sonic toothbrush:

Clean	82.000 mpm* – Efficiently removing dental plaque (2min)
White	96.000 mpm – Getting rid of surface stains and whitening teeth (2min)
Polish	90.000 mpm – Brightening and polishing front teeth (1min)
Face Care	70.000 mpm – Stimulating massage and cleansing the face skin (3min)

Sensitive – 62.000 mpm – An extremely gentle mode for sensitive teeth and gums (2min)

*Mpm: movements per minute

9


How to clean your SEYSSO sonic toothbrush

Please regularly clean your SEYSSO sonic toothbrush, wash the toothbrush handle and head. If necessary, please use soft cloth and mild soap to clean its surface. Please do not use sharp articles to clean it. After each use remove the brush head. Carefully clean the handle top and the pin to remove the toothpaste sediment.

Particular attention should be paid to the places of the seals: near the pin and around the button, where the accumulation of toothpaste favors the growth of microorganisms. Dry the handle with a soft cloth.


When cleaning the travel case and charging station, unplug them. Use a dry cloth to remove any dirt from them. Do not clean the case when the UV light bulb is hot. Do not clean the brush head, handle, charging station or travel case in the dishwasher. Do not clean the case and the charging station under the tap.

Please do not clean the brush head, handle, charging base or travel case in dish-washing machine.

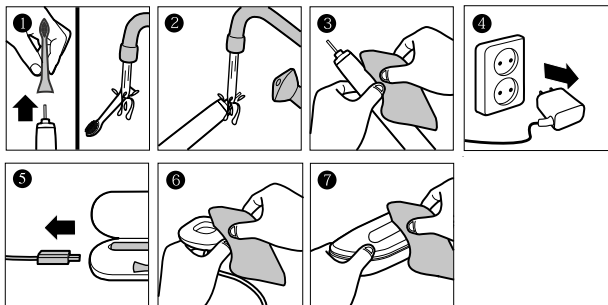10 How to store your SEYSSO sonic toothbrush

If the product is not to be used for a long time, unplug the charging station and travel case, clean all the components of the kit and store them in a cool, dry place away from sunlight.

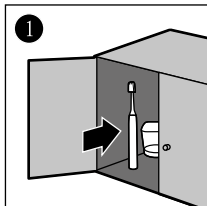11 Environmental protection

If you want to discard this product, please do not place it together with domestic waste but send it to the officially designated recycling center and dispose the battery at the officially designed battery recycling center. Doing so can benefit environmental protection.

The UV light bulb contains mercury that may be harmful to the environment. Do not break it and do not try to disassemble it.

12 Replacement heads

SEYSSO warranty will be invalid if it is found that the handle is damaged because non-SEYSSO replacement brush heads are used.

It is not recommended to use non-SEYSSO replacement brush heads, because SEYSSO cannot control the quality of other replacement brush heads.

SEYSSO cannot guarantee that non-SEYSSO replacement brush heads fit the handle.

Similarly to manual toothbrushes, sonic toothbrush heads should be replaced every 2-3 months.

The tips should be placed on the handle shaft at a right angle without using excessive force. The tip should remain a few millimeters from the brush. This is a deliberate solution that prevents two elements from hitting each other when vibrating.

The tips are a consumable element. As a result of pressing them against the teeth while brushing, minimal deformations of the plastic may occur, which will result in louder operation.

The difference in volume does not affect the operation of the handle. It should also not be audible when brushing (under gentle pressure) and after removing the tip.

All SEYSSO replacement brush heads are printed with the logo of SEYSSO and meet the high-quality standard of SEYSSO.

13 Warranty and services

The manufacturer provides a warranty for a period of two years from the date of purchase. The warranty is valid in Poland and on the premises of official distributors. The warranty covers damage resulting from manufacturing defects or incorrect workmanship of the product. Keep your proof of purchase for the warranty period. During the warranty period, you can submit a complaint if the product was used correctly and as intended.

THE WARRANTY DOES NOT COVER:

- 1) Surface damage resulting from normal use.
- 2) Damage resulting from:
 - a) accidental drop or impact, misuse, abuse or commercial use, negligence, modification of the device,
 - b) improper maintenance, carelessness of the product,

- c) connecting to the wrong power source,
- d) repairs or attempted repairs made by a party not authorized by the manufacturer. Any attempt to open or disassemble the device, its parts or accessories will void the warranty.

Warranty claims must be supported by evidence. The claim must be made within the warranty period. The material covering the handle is delicate, so be careful when using it. Scratches or abrasions on the surface of the handle will be treated as superficial damage resulting from normal use. Therefore, they will not be covered by the warranty. The warranty applies to components that affect the functioning of the inhaler. Proper care of your toothbrush significantly increases its lifespan. Failure to follow the recommendations may result in the loss of the warranty.

Input: 100-240 V~, 50/60 Hz

Power: 3 W

Rated Voltage: 4V/DC

This manual applies to model with numbers: SGFT10

PRODUCENT:

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WM Brands sp. z o. o.

www.seysso.pl

ul. Dziupli 12

02-454 Warszawa, Polska

Made in P.R.C.

The product complies with the requirements of the so-called "New Approach" Directives of the European Union (EU).

The manual and packaging are suitable for recycling.

This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU with the symbol of a crossed-out waste bin. This marking indicates that this equipment, after its service life, cannot be disposed of with other household waste. The user is obliged to hand it over to those conducting the collection of waste electrical and electronic equipment. Those conducting the collection including local collection points, stores, and municipal units, create a proper system that allows the disposal of this equipment. Proper disposal of waste electrical and electronic equipment helps avoid harmful consequences to human health and the natural environment, resulting from the presence of hazardous substances and improper storage and processing of such equipment.

Symbol indicating Class II protection for electrical devices.

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