**Oral-B**

**iO**

English	6
Polski	16
Český	28
Slovenský	39
Magyar	50
Hrvatski	64
Slovenski	75
Lietuvių	86
Latviski	97
Eesti	108
Українська	118
עברית	131

**Internet:**

[www.oralb.com](http://www.oralb.com)

[www.service.oralb.com](http://www.service.oralb.com)

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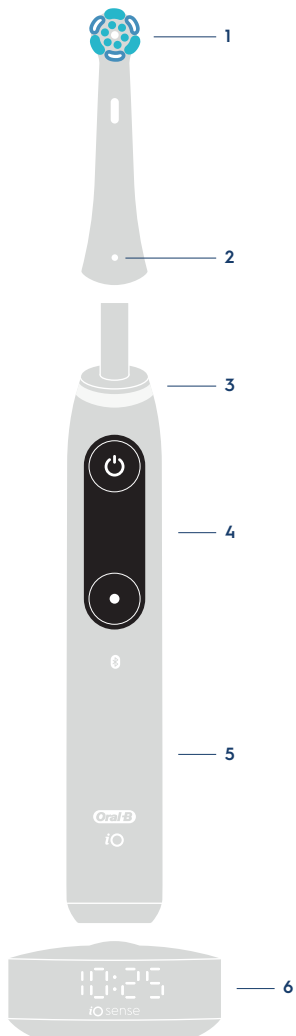
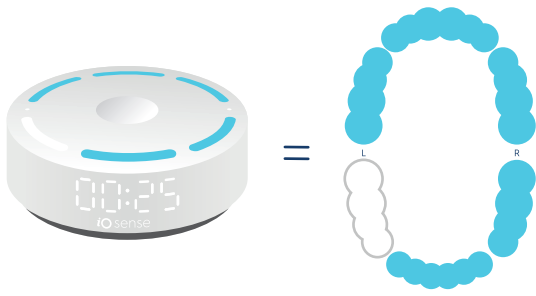
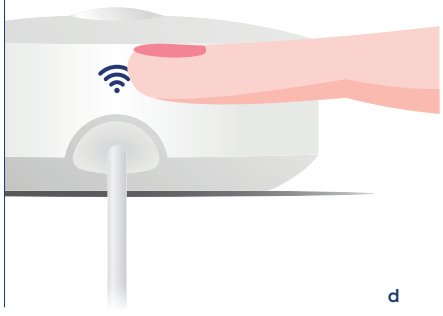
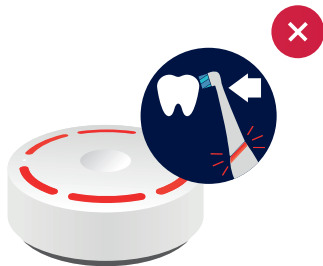
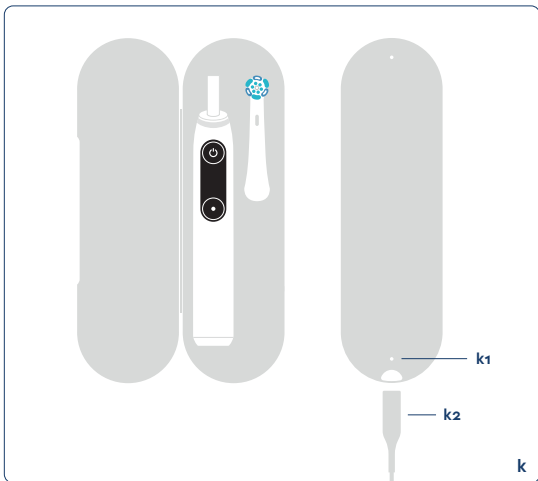
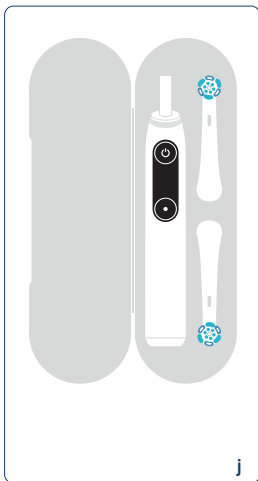
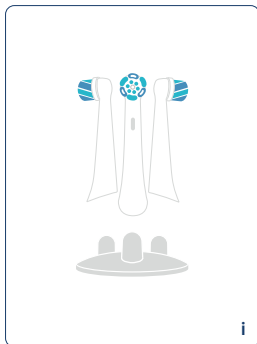
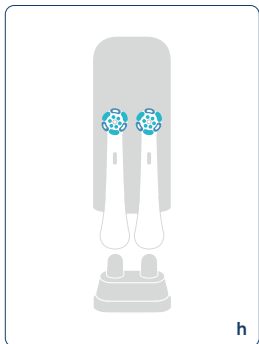
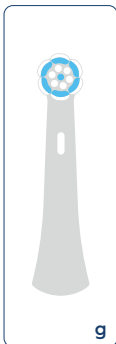
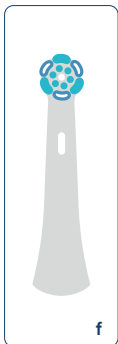
## Accessories (depending on model) and features:

- a Brushing recognition - Indicator LEDs
- b Brushing recognition - Brushing progress indication
- c Download and launch Oral-B App to enable all features
- d iOsense™ button
- e Smart Pressure Sensor
- f Brush head «Ultimate Clean»
- g Brush head «Gentle Care»
- h Brush head holder with cap
- i Brush head holder stand
- j Travel case
- k Charging travel case
- k1 Charging light
- k2 Power supply

Note: Content may vary based on the model purchased.

# English

Welcome to Oral-B!

Before operating your Oral-B iO toothbrush and iOsense™ charger, please read these instructions and save this manual for future reference.

## IMPORTANT

- Periodically check the entire product/cord/accessories for damage. A damaged or non-functioning unit should no longer be used. If the product/accessories is damaged, take it to an Oral-B Service Centre. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Never insert any object into any opening of the appliance.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without

adult supervision.

- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Only use the charger provided with your appliance.
- Only use the charger and/or special cord set with integrated Safety Extra Low Voltage power supply provided with your appliance. Do not exchange or tamper with any part of it, otherwise there is risk of an electric shock. If charging stand/charging travel case is marked with 492-XXXX, only use Braun/Oral-B special cord set with this identification.

## WARNING

- Do not place the charger, power supply or charging travel case in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Avoid sharp bending, jamming or squeezing of the cord.
- This appliance contains batteries that are non-replaceable. Do not open and disassemble device. For battery recycling please dispose whole device according to local environmental regulations. Opening the handle will destroy the appliance and invalidate the warranty.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is

not intended for use on multiple patients in a dental practice or institution. For hygienic reasons we recommend the use of individual brush handles and brush heads for each family member.

- Small parts may come off, keep out of reach of children.

To avoid brush head breakage which may generate a small parts choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fits properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Replace brush head every 3 months or sooner if brush head becomes worn.
- Clean brush head properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

## Important Information

- This toothbrush, charger and pouch (accessory) contain NdFeB magnets and generate magnetic fields. If you use a medical device or have an implanted medical device, consult with the medical device manufacturer or a physician, if you have questions about safely using your toothbrush.
- To avoid electromagnetic interference and/or compatibility conflicts, disable *Bluetooth* on your toothbrush handle (5) and *Bluetooth*/Wi-Fi on your iOsense™ charger before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals.
- Disable *Bluetooth* by entering the «Settings» menu on the interactive display (4) -> *Bluetooth* -> disable. Follow the same procedure to enable *Bluetooth* again. Disable *Bluetooth*/Wi-Fi on your iOsense™ charger by entering low power mode by tapping the iOsense™ button twice (d). Tap the button once again to enable *Bluetooth*/Wi-Fi.
- Persons with pacemakers should always keep the toothbrush or iOsense™ charger more than 15 centimeters (six inches) from the pacemaker while turned on. Anytime you suspect interference is taking place, disable *Bluetooth* on toothbrush and disable *Bluetooth* and *Wifi* on iOsense™ charger.

## Specifications

For voltage specifications refer to the bottom of the charging unit.

## First use

Please wake up your iO toothbrush by pressing the on/off button (7) and follow the instructions on the interactive display (4) to get started for first time setup.

If not further instructed on the display, you may toggle with the mode button (8) and select with the on/off button (7).

You may change the settings via the «Settings» menu (17) at any time.

You may connect your iO toothbrush with your iOsense™ charger by switch on the toothbrush in close proximity to the plugged-in iOsense™ charger.

## Charging and Operating

- To charge the toothbrush battery place the toothbrush on the iOsense™ charger (6).
- The charge level is indicated on the interactive display (4). The light ring (3) fades in and out during charging. Once it is fully charged the light ring turns off. A full charge typically takes 3 hours.
- If the battery is running low the light ring (3) blinks quickly in red. Once the battery is empty, the motor stops; it will need at least a 5-minute charge for one usage.  
Note: In case the battery is empty (display off), charge for at least 30 minutes.
- You can always store the toothbrush on the iOsense™ charger to maintain it at full charge; a battery overcharge is prevented by the device.

Note: Store the toothbrush at room temperature for optimal battery maintenance.

**Caution:** Do not expose the toothbrush to temperatures higher than 50 °C / 122 °F.

## Using your iO toothbrush

### Start Brushing

Place the brush head on the toothbrush handle. A small gap should remain between handle and brush head. Activate the display by picking-up the toothbrush, lifting it off the iOsense™ charger or pressing the on/off button (7) (display turns on). Wet brush head and apply toothpaste. To avoid splashing, guide the brush head to your teeth before switching the toothbrush on by pressing the on/off button.

Note: Display turns off after some time. You may press the on/off button to activate display again.

Guide the brush head from tooth to tooth. Brush

the outsides, the insides and the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you. After brushing, the display gives you feedback based on your brushing time.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» or «Super Sensitive» modes (optional in combination with an Oral-B «Gentle Care» brush head (g)). Availability of modes are depending on model purchased.

### **Toothbrush Timer (Stutter Timer)**

The iO toothbrush is equipped with timer features, which help you to achieve the professionally recommended 2-minute brushing time. The timer signals with a short stuttering sound at 30-seconds intervals to move to the next quadrant of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

The timer memorizes the reached brushing time, even when the toothbrush is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or placing the toothbrush on the plugged-in iOsense™ charger.

Note: During the use of the Oral-B App or iOsense™ charger features the stutter timer is deactivated in the toothbrush. You may also customize the stutter timer settings via Oral-B App.

### **Toothbrush Light Ring (3)**

Your iO toothbrush is equipped with a light ring (3) that has various functions indicated via colours:

#### **Smart Pressure Sensor (e)**

White (default) = Operating indicator (choose individual colour via Oral-B App) / Low pressure indicator

Green = recommended pressure indicator

Red = high pressure indicator

If the Smart Pressure Sensor illuminates in white (or your pre-selected colour), we recommend applying more pressure. The green light indicates recommended brushing pressure. If too much pressure is applied, it lights up in red reminding you to reduce pressure.

During high pressure the brush head movement will be reduced (in «Daily Clean», «Intense» and

«Sensitive» modes). Periodically check the operation of the Smart Pressure Sensor by pressing moderately on the brush head during use.

Note: During the use of the Oral-B App, the Smart Pressure Sensor indication may be inactive for some features.

Note: Smart Pressure Sensor is deactivated during «Tongue Clean» mode.

### **Other functions**

Red = Low battery indicator

Blue = Bluetooth pairing indicator

Orange = Error indicator (follow display messages)

### **Brushing modes (depending on model)**

Your iO toothbrush offers different brushing modes indicated on the display:

«Daily Clean» – Standard mode for daily cleaning

«Sensitive» – Gentle, yet thorough cleaning for sensitive areas

«Gum Care» – Gentle massage of gums

«Whiten» – Polishing for occasional or everyday use

«Intense» – For an extraordinary clean feeling (operated with a higher speed)

«Super Sensitive» – Super gentle cleaning for extra sensitive areas

«Tongue Clean» – Tongue cleaning for occasional or everyday use

Note: The sequence and availability of brushing modes may vary depending on model purchased.

When operating the «Tongue Clean» mode, we recommend using a «Gentle Care» brush head (g). You can brush your tongue either with or without toothpaste. Brush the whole area of your tongue systematically in gentle movements. The recommended brushing time is 30 seconds; the motor stops after 30 seconds.

#### **Mode setting:**

Wake up your iO toothbrush to identify selected mode for your brushing session. To change to others modes, successively press the mode button (8) during use. The toothbrush automatically starts in the mode selected last.

**Note:** You may also customise mode settings via Oral-B App.

### **Interactive Display**

Activate the interactive display (4) by picking-up the toothbrush, lifting it off the charger or pressing the on/off button (7).

(10-16) «Mode selection» – toggle by pressing the mode button during use.

(17) «Settings» – toggle through all modes until you reach the option «Settings»; to enter press the on/off button:

- (19) manage *Bluetooth* settings (enable, disable)
- (20) select/change colour of light ring
- (21) select/change language
- (22) perform factory reset
- (24) «Timer» – displays brushing time.
- (25) «Battery status» – shows battery status after use, when charging or when battery is low.
- (26) «Brush head replacement reminder» – reminds you when to replace your brush head.
- (27) «Brushing feedback» – brushing feedback is indicated with a smiley face.

When the toothbrush is picked up, lifted off the iOsense™ charger or when you press the on/off button the toothbrush handle indicates its readiness.

«Error messages» – «outside operating limits»: you are operating or charging outside standard temperature range.

## Brush heads

Oral-B offers you a variety of different Oral-B iO brush heads (f/g) that only fit your Oral-B iO toothbrush handle and can be used for tooth-by-tooth precision cleaning.

Most iO Oral-B brush heads feature INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the bristle colour will fade half-way approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the colour recedes, you may exert too much pressure on teeth and gums. In general, we recommend replacing the brush head every 3 months.

### Brush Head Replacement Reminder

The toothbrush is equipped with a «Brush Head Replacement Reminder» on the interactive display, indicating the need to exchange the brush head approximately after 3 months when used twice per day.

To reset the «Brush Head Replacement Reminder» follow the instructions on the display.

Note: You may also customize the brush replacement settings via the Oral-B App.

## Connecting your iO toothbrush to your Smartphone/Oral-B App

The Oral-B App offers various features:

- Personalize your toothbrush settings (e.g. choose your preferred light ring colour, adjust stutter timer, brush mode sequence or brush replacement reminder settings)
- Track your brushing behavior (e.g. brushing time and applied pressure)
- Real-time brushing recognition to help you brush every area for the right amount of time.

The Oral-B App is available for smartphones running iOS or Android™ operating systems. Check [app.oralb.com](http://app.oralb.com) for compatible smartphones and country availability of the Oral-B App. It can be downloaded for free from App Store<sup>(SM)</sup>/ Google Play store™ or you may scan the QR code on the cover page of this manual (c).

- Launch the Oral-B App. It will guide you through the entire *Bluetooth* pairing procedure. You may also pair new devices at any time in the Oral-B App's settings.

Note: Oral-B App functionality is limited if *Bluetooth* is deactivated on your smartphone (for directions please refer to your smartphone's user manual).

- You may change *Bluetooth* settings on your toothbrush via the display -> Settings -> *Bluetooth* enable/disable.
- Keep your smartphone in close proximity (within 5m distance) when using it with the toothbrush. Make sure the smartphone is safely positioned in a dry place.

Note: Your smartphone must support *Bluetooth* 4.2 (or higher) / *Bluetooth* Smart in order to pair your toothbrush.

**Caution:** Follow your smartphone's user manual to ensure your phone/ charger is designed for use in a bathroom environment.

## Using your iOsense™ charger

### Personalized live-coaching and Smart Pressure Sensor features

The iOsense™ charger has a built-in real-time feedback on brushing coverage, brushing time and pressure. The LED lights on top of your iOsense™ charger illustrate your upper and lower jaw (a). You may start your brushing session at any area of your mouth. While the area is being

brushed, the respective light fades from blue to white, indicating the brushing progress. Once one area is cleaned (a), the respective LED turns into white.

If too much pressure is applied, LEDs light up in red reminding you to reduce pressure (e).

**Caution:** During a brushing session with enabled brushing recognition feature on your iOsense™ charger, please stay within 5m distance to your iOsense™ charger to ensure a stable connection and brushing recognition.

To enable the full feature set of iOsense™ charger, the Oral-B App is required. Please follow the instructions in the Oral-B App for connection setup.

### Clock setup

Launch the Oral-B App and connect the iOsense™ charger with your Wi-Fi by following the instructions in the Oral-B App.

**Note:** The clock doesn't work without the Wi-Fi connection.

### Bluetooth & Wi-Fi® functions

The iOsense™ charger supports *Bluetooth* and *Wi-Fi*.

- *Bluetooth* connects the iOsense™ charger with the toothbrush and the Oral-B App.
- *Wi-Fi* connects the iOsense™ charger with the internet.


### iOsense™ button functions

You may find the capacitive touch button (d) on the back of your iOsense™ charger. It has two functions:

- **Low Power Mode:** Device enters into low power mode. All features beyond charging are disabled. Tap the button (d) twice to enable low power mode and once to disable. The activation of low power mode is indicated by the LEDs on the top of your iOsense charger. Power consumption of networked standby mode: 1,2 W.
- **Factory reset:** Touch and hold the button (d) until orange light appears on the top of your iOsense™ charger, then release the button. Tap the button once again while orange light is on, to trigger the factory reset. LEDs spin orange clockwise to confirm the factory reset. All credentials and settings are deleted from the iOsense™ charger.

## Cleaning Guidelines

After each brushing, rinse your brush head under running water with the handle switched on. Remove the brush head and clean the handle and brush head separately. Make sure water runs inside the top of the handle while the handle is switched on, as well as inside the brush head. Wipe all parts dry and let them air-dry before reassembling your toothbrush again. Unplug the charger before cleaning. Clean the charger regularly with a damp cloth. Please check individual cleaning instructions for each accessory in its «Accessories» paragraph.
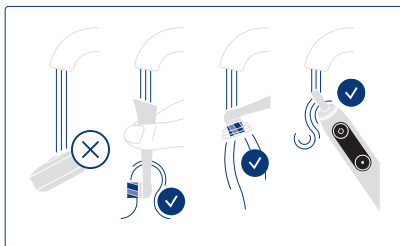
### **Accessories (depending on model)**

#### Charging Travel Case (k)

- To charge your toothbrush (5) place it in the charging travel case (k) and connect it to the mains outlet, using the extra power supply (k2) provided. The power supply may be used for all global mains voltage ranges (100 – 240V).
- The flashing light (k1) of the charging travel case indicates that the toothbrush is charging. Once the toothbrush is fully charged the light turns off. A full charge in the charging travel case typically takes 6 hours.
- The outer and inner case should be cleaned with a damp cloth only. Make sure that the toothbrush and brush head are completely dry before storing in the travel case.

**Note:** Clean stains on travel case surfaces immediately. Store travel case in a clean and dry place.

#### Travel Case (j)

Store the clean and dry iO toothbrush/brush heads in the travel case (j). The outer and inner case should be cleaned with a damp cloth only.

**Note:** Clean stains on travel case surfaces immediately. Store travel case in a clean and dry place.

### Brush head holder stand (i) / Brush head holder with cap (h)

You may store your cleaned and dry brush heads on the pins of the brush head holders. Make sure your brush heads are dried before placing them in the holder and/or closing the cap. Brush head holders and cap are dishwasher-safe.

### Zipper Pouch / Magnetic Pouch

You may store your cleaned brush heads and toothbrush in the pouches. Make sure your brush heads and toothbrush are dried before placing them in the pouches. The outer and inner surfaces can be cleaned under running water. Store pouches in a clean and dry place.

Note: Clean stains on pouch surfaces immediately.

Subject to change without notice.

## Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

## Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. Depending on availability replacement of a unit might result in a different colour or an equivalent model.

This guarantee extends to every country where this appliance is supplied by Braun / Oral-B or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun / Oral-B parts are not used.

These appliances are equipped with a *Bluetooth* Smart approved radio/Wi-Fi interfaces. Failure to establish a *Bluetooth* Smart/Wi-Fi connection to

specific smartphones is not covered under the device guarantee, unless the malfunction of the *Bluetooth*/Wi-Fi interface is referable to Oral-B. The *Bluetooth*/Wi-Fi devices are warranted by their own manufacturers and not by Oral-B. Oral-B does not influence or make any recommendations to device manufacturers, and therefore Oral-B does not assume any responsibility for the number of compatible devices with our *Bluetooth*/Wi-Fi interface.

Oral-B reserves the right, without any prior notice, to make any technical modifications or changes to system implementation of device features, interface, and menu structure changes that are deemed necessary to ensure that Oral-B systems function reliably.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

### Information

#### ***Bluetooth*® radio module/ Wi-Fi® module**

Although all indicated functions on the *Bluetooth*/Wi-Fi device are supported, Oral-B does not ensure 100% reliability in connection, and feature operation consistency.

Operation performance and connection reliability are direct consequences of each individual *Bluetooth*/Wi-Fi device, software version, as well as the operating system of those *Bluetooth*/Wi-Fi devices and company security regulations implemented on the device.

Oral-B adheres to and strictly implements the *Bluetooth*/Wi-Fi standard by which *Bluetooth*/Wi-Fi devices can communicate and function in Oral-B iO toothbrushes.

However, if device manufacturers fail to implement this standard, *Bluetooth*/Wi-Fi compatibility and features will be further compromised and the user may experience functional and feature-related issues. Please note that the software on the *Bluetooth*/Wi-Fi device can significantly affect compatibility and operation.

### Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable toothbrush handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.

- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

## Troubleshooting

Problem	Possible reason	Remedy
<b>Oral-B App</b>		
Toothbrush and iOsense™ charger do not operate (properly) with Oral-B App.	<ol style="list-style-type: none"> <li>1. Oral-B App is closed.</li> <li>2. <i>Bluetooth</i> is not activated on smartphone.</li> <li>3. <i>Bluetooth</i> is disabled on toothbrush</li> <li>4. <i>Bluetooth</i> connection to smartphone has been lost.</li> <li>5. Your smartphone does not support <i>Bluetooth</i> 4.2 (or higher) / <i>Bluetooth</i> Smart.</li> <li>6. Outdated Oral-B App.</li> <li>7. Your smartphone is not compatible with the Oral-B App.</li> <li>8. iOsense™ charger is in Low Power mode.</li> </ol>	<ol style="list-style-type: none"> <li>1. Launch Oral-B App.</li> <li>2. Activate <i>Bluetooth</i> on smartphone (described in its user manual).</li> <li>3. Enable <i>Bluetooth</i> via display settings on your toothbrush.</li> <li>4. Keep your smartphone in close proximity when using it with the toothbrush and iOsense™ charger.</li> <li>5. Your smartphone must support <i>Bluetooth</i> 4.2 (or higher)/ <i>Bluetooth</i> Smart in order to pair toothbrush.</li> <li>6. Download latest Oral-B App.</li> <li>7. Check <a href="http://app.oralb.com">app.oralb.com</a> for compatible smartphones and country availability of the Oral-B App.</li> <li>8. Quit Low Power mode on iOsense™ charger by tapping the button (d) once.</li> </ol>
<b>TOOTHBRUSH</b>		
Unwanted toothbrush behavior	Toothbrush was misconfigured	Perform a factory reset to restore initial settings. Navigate on toothbrush display to Settings -> Factory Reset (22). Press on/off button to confirm.
Toothbrush is not operating (during first use).	<ol style="list-style-type: none"> <li>1. Battery has a very low charge level; handle does not operate.</li> <li>2. Operating outside standard temperature range (message on display and/or alternating lights of light ring)</li> </ol>	<ol style="list-style-type: none"> <li>1. Charge at least for 30 min.</li> <li>2. Always operate in standard operating temperature range (&gt; 0 °C and &lt; 40 °C).</li> </ol>

Short stuttering sounds after 2 minutes or toothbrush does not provide any stutter timer.	<ol style="list-style-type: none"> <li>1. Timer has been modified/ deactivated via Oral-B App.</li> <li>2. The toothbrush timer is disabled during use of brushing recognition with Oral-B App or iOsense™ charger.</li> </ol>	<ol style="list-style-type: none"> <li>1. Use Oral-B App to change timer settings or perform toothbrush factory reset.</li> <li>2. Disconnect toothbrush from Oral-B App and iOsense™ charger or configure iOsense™ charger with the App to disable brushing recognition.</li> </ol>
Toothbrush does not charge.	<ol style="list-style-type: none"> <li>1. Battery is at a very low charger level (0%); it may take up to 10-15min until display turns on.</li> <li>2. Ambient temperature for charging is out of valid temperature range (&lt; 0°C and &gt; 40°C).</li> <li>3. Charger may not be (properly) connected.</li> </ol>	<ol style="list-style-type: none"> <li>1. Charge at least for 30 min.</li> <li>2. Recommended ambient temperature range for charging is 5°C to 35°C.</li> <li>3. Ensure iOsense™ charger is connected to mains voltage.</li> </ol>
Rattling/loud noise during brushing or turning on.	<ol style="list-style-type: none"> <li>1. Brush head may not be placed correctly onto the brush handle.</li> <li>2. Brush head may be worn out.</li> </ol>	<ol style="list-style-type: none"> <li>1. Check proper fit of brush head on handle. Small gap remains in between.</li> <li>2. Replace brush head.</li> </ol>
<b>iOsense™ charger</b>		
Clock is not displayed on the iOsense™ charger or an animation is playing ongoingly on clock display	Device is not connected to the internet.	Use the Oral-B App to setup the iOsense™ charger. Please follow the instruction in the Oral-B App. Make sure your Wi-Fi network is in range of iOsense™ charger and connected to internet.
iOsense™ charger is not reacting when using my toothbrush or reacts on a different toothbrush	<ol style="list-style-type: none"> <li>1. iOsense™ charger is connected to a wrong toothbrush</li> <li>2. iOsense™ charger is in Low Power mode.</li> </ol>	<ol style="list-style-type: none"> <li>1.1 When you are not using the Oral-B App: <ul style="list-style-type: none"> <li>- Enable <i>Bluetooth</i> via display settings on your toothbrush</li> <li>- Reset iOsense™ charger.</li> <li>- Hold your toothbrush close to iOsense™ charger and turn it on to connect both devices.</li> </ul> </li> <li>1.2 When you are using the Oral-B App: <ul style="list-style-type: none"> <li>- Enable <i>Bluetooth</i> via display settings on your toothbrush.</li> <li>- Use the Oral-B App to connect your toothbrush.</li> </ul> </li> <li>2. Quit Low Power mode by tapping the iOsense™ button (d) once.</li> </ol>
iOsense™ charger is not showing my current brushing session	<ol style="list-style-type: none"> <li>1. iOsense™ charger may got disconnected from your toothbrush while brushing</li> <li>2. iOsense™ charger is in Low Power mode.</li> <li>3. Brushing feedback is already provided in the Oral-B App</li> </ol>	<ol style="list-style-type: none"> <li>1. Please ensure that you stay in close proximity to your charger while brushing</li> <li>2. Quit Low Power mode by tapping the iOsense™ button (d) once.</li> <li>3. Brushing feedback can only be provided on one device. Make sure you have selected the right device.</li> </ol>

The rechargeable iO toothbrush contains a radio module which operates in the 2.4 to 2.48 GHz band with a maximum power of 1mW.

The iOsense™ charger contains a radio module which operates in the 2.4 to 2.48 GHz band with a maximum power of 100mW.

Hereby, Braun GmbH declares that the radio equipment types 3786, 3768, 3758, 3771, 3765, 3754, 3764, 3762 and 3767 are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address:  
[www.oralb.com/ce](http://www.oralb.com/ce).

#### Manufacturing date code:

To determine the date of manufacture, refer to the production code located on the type plate on the bottom of the **toothbrush handle** and/or **iOsense™ charger**.

#### Toothbrush Handle:

The production code starts with two letters. The first digit after the two letters of the production code refers to the year of manufacture. Second and third digits refer to the calendar week of manufacture. Fourth digit refers to the calendar day of manufacture.

Handle: Example: «XX1093...» – means the product was manufactured on Wednesday in week 09 of 2021.

#### iOsense™ charger:

The production code of the device consists of 5 digits. The first three digits represent the Julian Date (Example: 035 = 4 February), the fourth digit is year of production (Example: 1 = 2021), the fifth digit is an internal code.

#### Charging travel case:

The production code of the device consists of 6 digits. The first two digits represent the Year (Example: 21 = 2021), the next 2 digits are the calendar week of production (Example: 01-52), the last two digits are internal codes.

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