
# AQUACARE WATERFLOSSER

powered by **BRAUN**

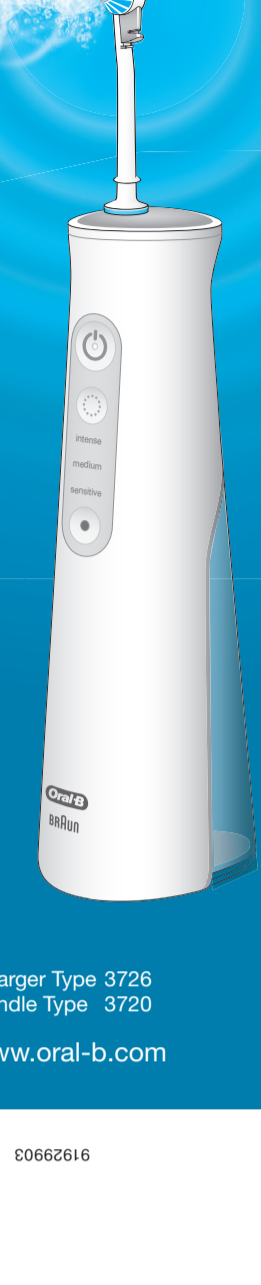Charger Type 3726  
Handle Type 3720

[www.oral-b.com](http://www.oral-b.com)

80662616

**UK** 0 800 731 1792

**IE** 1 800 509 448

**PL** 801 127 286  
801 1 BRAUN

**CZ** 800 11 33 22

**SK** 800 333 233

**HU** (06-1) 451-1256

**HR** 01 / 66 90 330

**SI** 080 2822

**LT** (8 5) 205 1272

**LV** 67798667












**EE** 667 5047

**Internet:**  
[www.oralb.com](http://www.oralb.com)  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)

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91929903/X-18  
EN/PL/CZ/SK/HU/HR/SI/LT/LV/EE
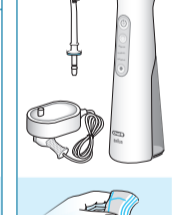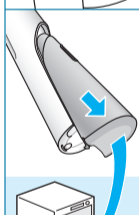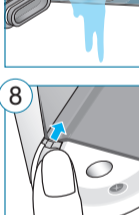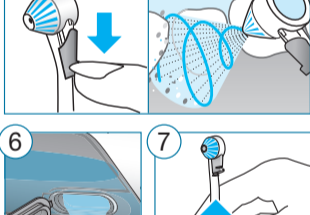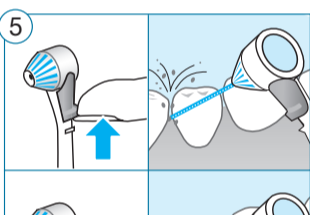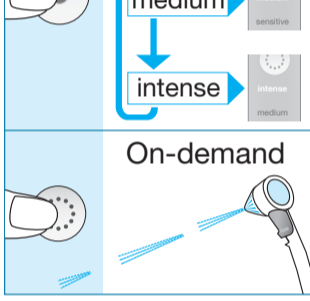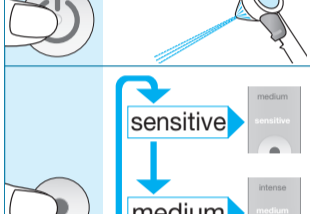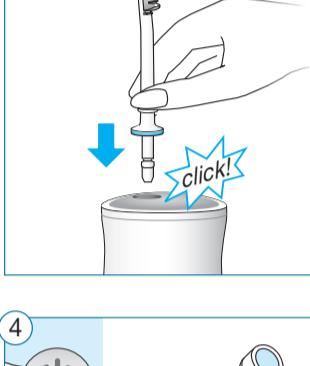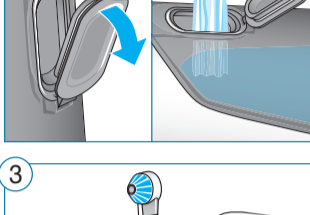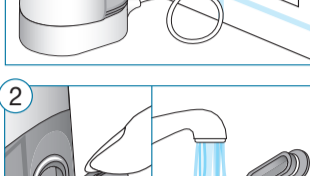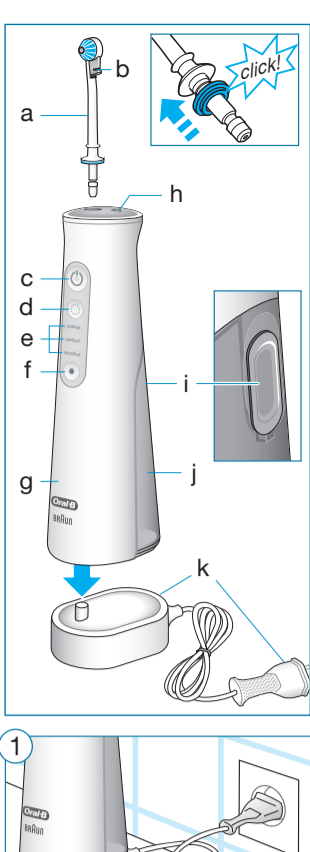
## English

Welcome to Oral-B!  
Please read the usage instructions carefully before first use and save this manual for future reference.

### WARNING

- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable. Do not disassemble the product except when disposing of the battery. When taking out the battery from the handle for disposal, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This appliance is a personal care device and is not intended for use on multiple persons.
- Small parts may come off, keep out of reach of children.

## Using the irrigator

Fill the water reservoir (j) with lukewarm water by lifting the reservoir cap (i) on the back of the handle (g). Lean over the sink and direct the nozzle in your mouth. Close your lips slightly and turn on the handle by pressing the on/off button (c). Guide the jet along the gumline and pause briefly between teeth. Alternatively, you may also point the nozzle tip at a desired area and use the on-demand button (d) for a single jet stream. Keep your mouth slightly open to allow the water to run into the sink. Refill water reservoir as many times as needed to complete your cleaning routine. It may take multiple refills to clean your entire mouth.

## After use

Make sure the water reservoir is completely empty. Press the on/off button (c) and leave it turned on until no more water is running through the nozzle, to prevent any deposit build-up or bacteria from collecting in the standing water. Turn off the handle by pressing the on/off button again once the water is empty. A longer operation of the handle without water may damage the motor. To drain any water out of the handle hold it over the washbasin and tilt it from side to side (picture 6).

Release the nozzle (a) by pressing the release button (h) and pulling the nozzle at the same time (picture 7).

Lift up the cap (i) and let the water reservoir air-dry. Always wipe the handle and nozzle dry. Note: Make sure you always empty the water reservoir before storing or travelling with the irrigator.

In general we recommend to replace the nozzle every 3 months.

## Cleaning recommendations

Clean and dry handle after each use and leave the water reservoir cap (i) open. The water reservoir (j) is dishwasher-safe and can be detached via the snap-fit on the bottom of the handle. Move the lever and slide the reservoir down (picture 8). Regularly clean the charger and handle with a damp cloth. Store the appliance in a cool and dry place.

Subject to change without notice.

## Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

## Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads/nozzles as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

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